

NISTIR 88-3876



An Investigation of a User-Operated Mass Calibration Package

Randall M. Schoonover and James E. Taylor

U.S. DEPARTMENT OF COMMERCE
National Institute of Standards and Technology
(Formerly National Bureau of Standards)
National Engineering Laboratory
Center for Manufacturing Engineering
Gaithersburg, MD 20899

November 1988

Prepared for
U.S. DEPARTMENT OF COMMERCE
National Institute of Standards and Technology
(Formerly National Bureau of Standards)
Office of Weights and Measures
Gaithersburg, MD 20899



75 Years Stimulating America's Progress
1912-1988

NISTIR 88-3876

An Investigation of a User-Operated Mass Calibration Package

Randall M. Schoonover and James E. Taylor

U.S. DEPARTMENT OF COMMERCE
National Institute of Standards and Technology
(Formerly National Bureau of Standards)
National Engineering Laboratory
Center for Manufacturing Engineering
Gaithersburg, MD 20899

November 1988



National Bureau of Standards became the National Institute of Standards and Technology on August 23, 1988, when the Omnibus Trade and Competitiveness Act was signed. NIST retains all NBS functions. Its new programs will encourage improved use of technology by U.S. industry.

Prepared for
U.S. DEPARTMENT OF COMMERCE
National Institute of Standards and Technology
(Formerly National Bureau of Standards)
Office of Weights and Measures
Gaithersburg, MD 20899

U.S. DEPARTMENT OF COMMERCE
C. William Verity, Secretary
NATIONAL INSTITUTE OF STANDARDS
AND TECHNOLOGY
Ernest Ambler, Director

An Investigation of a User-Operated Mass Calibration Package

by

Randall M. Schoonover and James E. Taylor

Abstract

Reported here are the results of two round-robin mass measurement programs. The first round-robin elucidated the many technical problems that required solutions before successful mass calibrations could be performed in all of the participant laboratories. This report details the technical innovations incorporated in the second round-robin, i.e., thermal conditioning of the kilogram weights, balance servo control, automatic data acquisition, the measurement of some air density parameters, and computer software and presents the results. We believe these results clearly indicate that properly instructed personnel could, in the future, successfully calibrate mass standards at locations remote from the National Institute of Standards and Technology (NIST) laboratory while maintaining the rigor necessary for certification.

Introduction

The process of disseminating the mass unit begins when a metrologist observes the difference between two masses by noting the force each mass exerts on a balance beam. One of the masses is the standard and embodies the mass unit, and the other is to be assigned a value relative to the standard. In practice more information is required on both objects than simply the mass of the standard. The metrologist needs additional information that characterizes both objects and, in addition, detailed knowledge of the entire calibration process.

In essence the work reported here is the application of knowledge gained from planned studies of systematic errors in the mass measurement process. Here we describe a user-operated mass calibration package whose design builds on much of the early work. We also provide data from a ten laboratory round-robin test of the package. We believe the test results support the concept of a mass calibration performed by the user and certified by the NIST.

Furthermore, we believe the concept of this package is unusual in that it contains everything necessary to perform the calibration, insure its correctness and issue a calibration report. In short, all the user provides is laboratory space, a balance of appropriate capacity and precision, and a balance operator.

Philosophy

In considering the usefulness of a user-performed mass calibration, one could compare it to the established service where the user submits mass standards to the NIST for calibration. Other than a purchase order no

other communication is necessary between the user and the NIST. In short, the user pays a fee and the NIST provides a report of calibration. The approach we propose offers user participation in the calibration process, enhances communication and, according to most users leaves them with much more than a report of calibration. In summary, they gain confidence in the measurements they provide to their clientele as this method disseminates techniques and experience as well as the mass unit.

Background

Beginning in the early 1960's the NIST mass group became interested in the development of a solid object density scale. The thinking at that time was that the material would probably have a density lower than stainless steel and be either of single crystal in structure or a ceramic [1]. The material eventually chosen was single crystal-silicon with a density of 2.3 grams per cubic centimeter. The choice of silicon immediately posed many questions as our experience in high accuracy weighing was limited to plated brass or stainless steel laboratory mass standards. There was a concern regarding the adequacy of the air density equation used for making air buoyancy corrections. Such corrections would be unusually large when making mass comparisons between silicon and metallic weights with densities near 8 grams per cubic centimeter. The correction is almost non-existent when weighing various alloys of stainless steel. There were other concerns about weighing dissimilar materials related to surface effects, stability, thermal and electrical conductivity, etc. Interest in these matters was again renewed in 1975 by Pontius [2] and has continued until present. We have summarized all of this work in a paper published in 1986 [3]. However, the individual studies can be reviewed in detail by referring to references [2-8]. For completeness we give whatever brief descriptions are necessary as the need arises in the text of this report.

Round-Robin Measurements

In the report of our earlier round-robin measurement [7] we suggested that any future package should include the following:

- o Portable computer with data reduction software.
- o Provision for thermal conditioning of the masses.
- o Calibrated instruments for measuring the air density parameters.
- o Automation to reduce work load.
- o A pair of calibrated kilogram mass standards.

In the work reported here we have successfully provided for all of the above except for the one-kilogram mass standards. In addition we provided more automation than was anticipated in 1983 when the earlier round-robin was completed. Beyond these initial goals we also included more laboratories at greater distances to test the equipment fragility and to periodically check the correctness of the instrument calibrations.

The laboratories included in this round-robin in sequence are:

- 1) May 86 - The National Bureau of Standards-Gaithersburg, Maryland.
- 2) June 86 - State Laboratory-Sacramento, California.
- 3) July 86 - State Laboratory-Honolulu, Hawaii.
- 4) August 86 - 3M Laboratory-St. Paul, Minnesota.*
- 5) September 86 - State Laboratory-Reno, Nevada.
- 6) October 86 - State Laboratory-Salt Lake City, Utah.
- 7) November 86 - State Laboratory-Denver, Colorado.
- 8) December 86 - Sandia Laboratory-Albuquerque, New Mexico.*
- 9) May 87 - State Laboratory-Albany, New York.
- 10) June 87 - Mound Laboratory-Miamisburg, Ohio.*
- 11) July 87 - The National Bureau of Standards-Gaithersburg, Maryland.

*Laboratories capable of performing instrument calibration checks.

The measurements performed during this round-robin are identical to those of the earlier work with the addition of an extra kilogram pair. The five kilograms in the package are designated A1, S2, R2, H2, and TE. Although only one kilogram was added to the package, the specified weighing technique requires the kilograms be used in pairs; as a consequence kilogram A1 is used twice. The pairings were (A1,S2), (A1, TE), and (R2,H2). A brief recapitulation of the weighing scheme is given here for the convenience of the reader. Three series of measurements are performed in which each kilogram pair in the package is compared to the participating laboratory's pair of kilogram mass standards by means of substitution weighings on a single pan balance. In each series, six pairwise comparisons are made between the four kilograms involved. The resulting data are adjusted by the method of "least squares" and mass values are assigned to the kilogram pair in the package. Additionally this weighing technique yields statistical information about the weighing process and the standards. As there are three pairs of kilograms in the package, there must be three such groups of weighings performed by each participant. This method of intercomparing four like objects is referred to as a "four ones series" and is commonly used by the NIST mass calibration service.

Balance and Balance Servo

A. Balance

Each participating laboratory has in its possession a single-pan, one-kilogram capacity, substitution balance of adequate precision for routine mass calibration. All of these balances are made by the same manufacturer with some being more modern than others. However, all of them perform to the same specifications. Unfortunately, these balances are all mechanical and lack the electronic output signal necessary for automatic data acquisition. The balance weighing data are observed visually and then entered into the computer manually with the exception of the NIST balance that has been modified as described in the next section.

B. Balance Servo

At this time there is no general-purpose, single-pan balance of 1-kilogram capacity, with adequate precision for this work, and with an electronic data output. Robert D. Cutkosky, formerly of the Length and Mass Division, graciously consented to design an electronic balance servo that provided our balance with an electronic signal output. We assembled the circuit and modified the balance for use with the servo, with such success that we also installed one on the primary mass comparator here at NIST (the Voland Balance¹). Since this is an important new addition to this project and will also be retrofitted to other balances used in the mass calibration program we will describe it in detail.

Figure 1 is the schematic diagram of the balance servo circuit. We believe that with modification of the feedback components this servo can be installed on other similar balances but we have only demonstrated this on balances with a one kilogram capacity and arm lengths of approximately 5 cm. The balance was modified to permit the mounting of an electromagnet coil and photodetector mask to the balance beam as depicted in Fig. 2. Other modifications were performed that allowed for the mounting of permanent magnets, a photosensitive error detector using a pair of matched diodes, and electric conductors connected to the coil on the beam.

In principle the servo applies a magnetic force to the balance beam that maintains the relative position between the error detecting diodes rigidly mounted to the balance case and the photo mask attached to the free swinging beam. The force is generated by the circuit driving electrical current through the coil and the moving charges react with the fields of permanent magnets causing the desired motion of the balance beam. The circuit provides a current that is proportional to the displacement, velocity, and the time-integrated displacement of the beam whenever the photo mask is off its null position. The null position is selected by adjusting the free swinging rest point of the beam to be at the gravitational horizon. This minimizes coupling of pan swing to the beam [9] and also the offset current through the coil.

¹Certain trade names and company products are identified in order to adequately specify the experimental procedure. In no case does such identification imply recommendation or endorsement by the National Institute of Standards and Technology nor does it imply that the products are necessarily the best available for the purpose.

Figure 1. Schematic of the balance servo circuit.


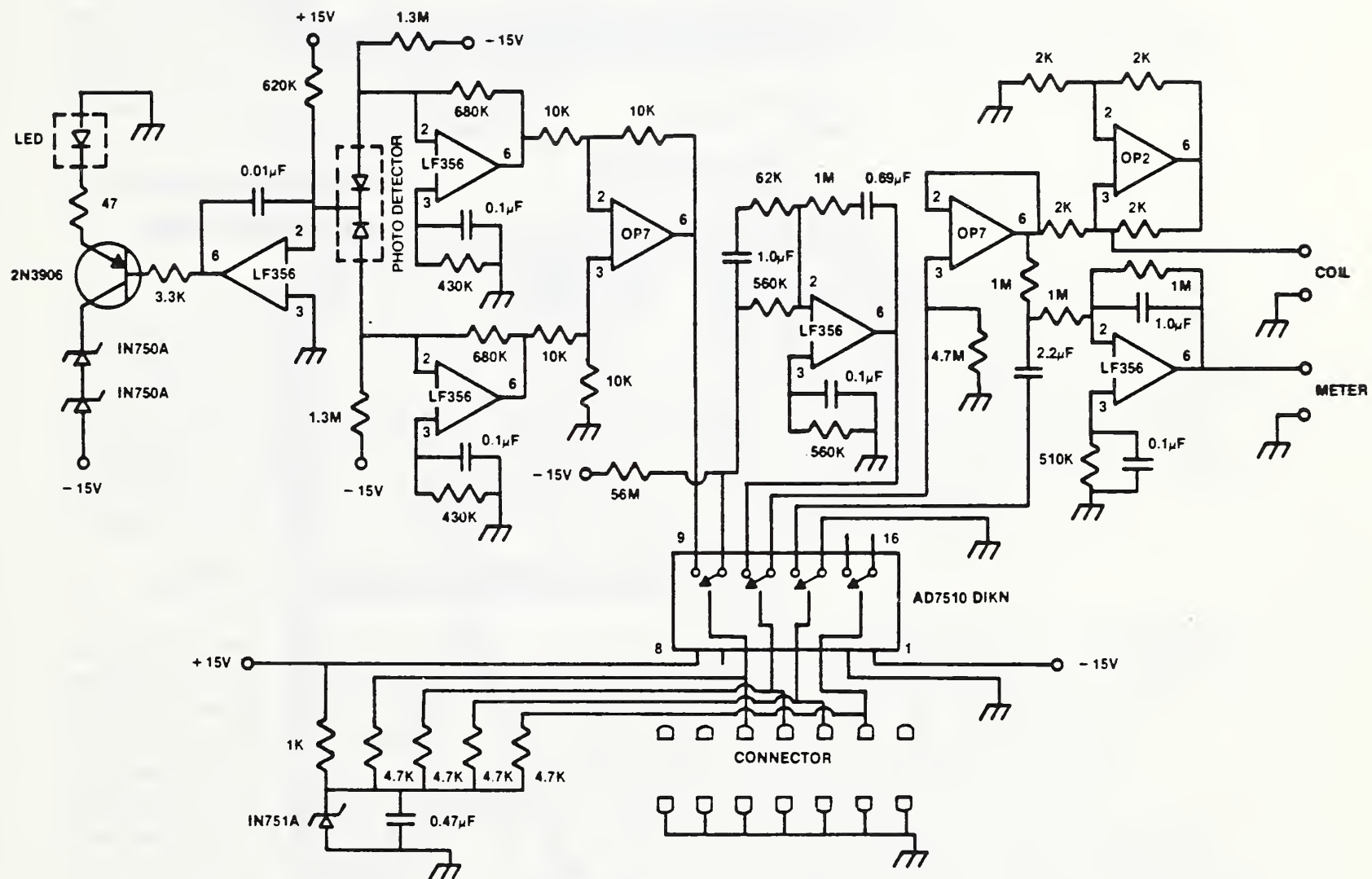


Figure 2. Configuration of balance and added servo components.

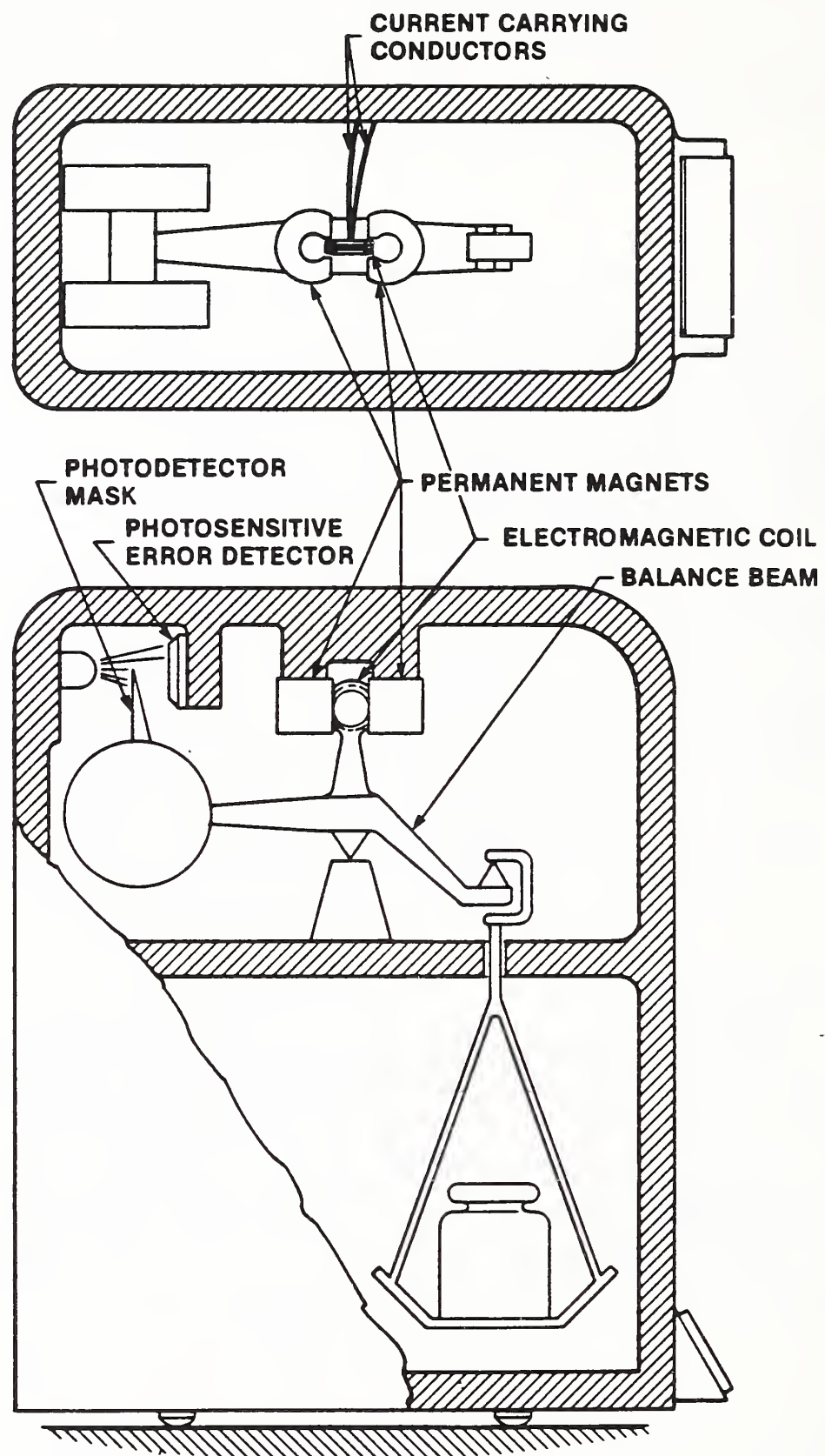

Initially the servo was tested with the balance modified to operate under constant load using a weight changing mechanism that we installed. Operating the balance in this configuration required other extensive modifications that are, in principle, described in [10]. The results indicate a significant improvement in precision, a constant sensitivity, and a larger on-scale range. Test data indicate we could tolerate differences four times larger than the original 100 milligrams and the balance response remained linear. We found it cumbersome to use the weight changing mechanism in this work and resorted to the much faster hand loading method. In doing so the precision reverts to that of the unservoed balance but the other virtues, especially the electronic data output that makes computer data acquisition possible, remain. Table 1 summarizes the precision tests. Although the balance reproducibility was not improved when operated in the normal manner, where the knife-flat contact is broken between each load, this may not be the case for every balance. Our balance precision had deteriorated slightly from its former level, perhaps to a degree and from a source that any improvement from the servo would not be evident.

Experimental Package

A. Mass Artifacts

The mass artifacts of the package are the same as those described in our earlier work with the exception of an addition of a one-kilogram weight designated TE. In brief, all of the masses are nominally one kilogram and are made of stainless steel. The densities are in the range of 7.8 to 8.0 g/cm³ except for weight H2 which is hollow and has a density of 2.9 g/cm³. Likewise, all of the weights are of one piece construction except for S2 which is comprised of three concentric rings and a supporting base. Table 2 summarizes the pertinent characteristics for each weight.

B. Thermally Controlled Soaking Plate.

Based on our study of surface-dependent thermal effects [8] it had become evident that systematic errors arise from the lack of thermal equilibrium between weights undergoing calibration and the surrounding air in the balance weighing chamber. To remedy this situation we modified a temperature control circuit designed by Cutkosky and Davis [11] to maintain the weight-storage area outside the balance at the same temperature found inside the balance weighing chamber. A heat source that has about the same effect on the balance as the operator's body heat is also part of this solution. In essence the temperature controller has two platinum thermometer sensors, one inside the balance and one buried in a nearby aluminum soaking plate where the weights are stored. Thus the soaking plate and the weights on it are in close proximity to the temperature in the balance weighing chamber. The rise in temperature associated with the operator's presence is anticipated and simulated by an electrically resistive heat source. This source does not radiate in the visible range and is replaced by the operator during the weighing period. The soaking plate is fabricated from solid aluminum alloy 2.54 cm thick and is 21.6 cm by 45.7 cm in plane view. The controlled heating element

Table 1. Comparison of balance precision using a balance in two different modes of operation.


Balance Configuration	Standard Deviation (μ g)		
	Test #1	Test #2	Test #3
o servoed beam o constant load mode	6	17	18
o automatic weight changer			
o unservoed beam o knife arrestment between loads	52	25	44
o hand loading of weights			

Table 2. Physical characteristics of each weight.

	Kilogram Designation				
	A1	S2	R2	H2	TE
Mass (g)	1000	1000	1000	1000	1000
Volume (cm ³)	127	126	126	338	125
Density (g/cm ³)	7.85	7.90	7.91	2.96	7.99
Surface Area (cm ²)	145	660	270	270	180
Height Ratio (H _x /H _{A1}) ¹	1.0	1.2	1.4	1.4	3.5

1. H_x (H_{A1}, H_{S2}, H_{R2}, H_{H2}, or H_{TE}) is defined as the height of each weight.

Figure 3. Soaking plate temperature control circuit.

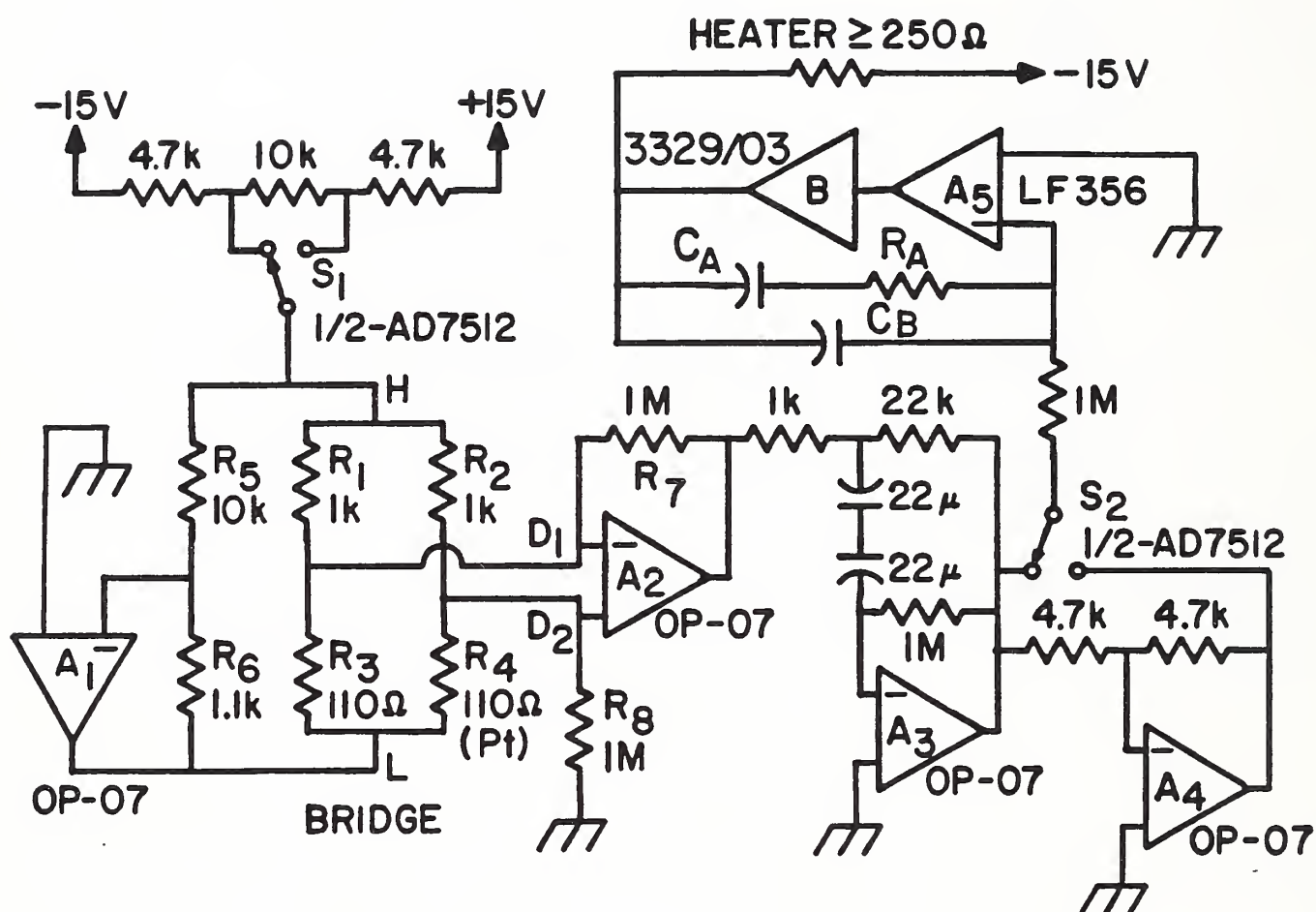

is a series hookup of 16 resistors of 16.2 ohms each. We make the following comments about the analog portion of the circuit shown in Fig. 3.

- 1) A 10k ohm resistor was added between switch S1 and the junction of resistors R1 and R5.
- 2) Resistor R3 is replaced with a 110-ohm platinum thermometer (R3 and R4 are both industrial quality thermometers).
- 3) The circuit is balanced by shunting the bridge with resistors, as necessary, to attain zero heater power when the platinum thermometers are in thermal equilibrium with each other. Shunt resistors less than 100k ohms are mounted on the aluminum plate to enhance the thermal stability of the circuit.
- 4) Capacitors C_A , C_B , and resistor R_A are respectively 5 μF , 0.022 μF , and 56 megohms. These components allow the controller to regulate the temperature of a soaking plate as described above including the thermal insulation.

C. Thermometer

A dual probe thermistor thermometer was selected for this work because of its small probe size, speed of response, and an optional analog electrical signal output that is needed for computer data acquisition. Unfortunately the probe selection can not be computer controlled on this instrument; the selection is made manually by the operator. By experiment it was determined that the air temperature, t , could be calculated from the analog voltage output using the following relationship:

$$t = 84.02531V - 0.101,$$

where V is the analog voltage output. For measuring air temperature one of the thermometer probes is suspended inside the balance case near the weighing pan and the other probe is buried in the aluminum soaking plate outside the balance where the weights are stored for thermal soaking before weighing.

The thermometer probe inside the balance can be monitored by the computer before and during the weighing procedure. This feature is important to assure that the thermal history prior to weighing is conducive to a successful weighing outcome. With calibration the thermometer is accurate to within 0.02 °C for either probe. Accuracy of the thermometer is achieved by calibration of both probes at several temperatures between 0 and 25 °C. This is accomplished by comparing the probes to a thermometer that is estimated to have an uncertainty in temperature of 0.003 °C. Field checks on the thermometer can be obtained by comparing the probes to a mercury-in-glass thermometer inside the balance weighing chamber, or by comparing the difference between the probes themselves when they are

thermally soaked at the same temperature; the observed difference can be compared to an established value, providing the user some assurance that the thermometer has retained its calibration.

D. Barometer

Of the three parameters that must be measured to compute air density -- temperature, pressure, and relative humidity -- in practice, pressure is the most difficult to measure. Many laboratories do not have the means to define pressure and must rely on instruments that require calibration. Without a local source of calibration, laboratories are forced to rely on rather delicate instruments to retain their calibration in transit from distant facilities.

The needs of this measurement project are well satisfied with an accuracy of a few parts in ten thousand. We selected a capacitance pressure transducer with an electronic signal output for ease in automatic data acquisition (0 to 10 volts) and with a pressure range from 600 mm Hg to 825 mm Hg. We anticipated a pressure range from 615 mm Hg to 760 mm Hg related to the range in elevation of the participating laboratories: from sea level to 3000 m. The relationship between the voltage output, V, and pressure, P, is as follows:

$$P = \frac{(30V+800)}{1.333224} + C$$

Where C is a small correction determined by calibration.

To test the ruggedness of the transducer before its inclusion in the package we shipped it to Sandia laboratories where it was calibrated; as it was before and after shipment at the NIST. With the normal precautions taken in the packing and shipping of delicate instruments we had no shift in the instrument calibration during the test or the round-robin measurements.

E. Relative Humidity Measurement

The measurement of relative humidity is accomplished with a thin polymer film capacitance probe that has a fast response time. Like the other instruments described above there is an analog signal output that is useful for automatic data acquisition. The relative humidity, RH, as a function of the voltage output is

$$RH = 1000V+C$$

where V is the voltage and C is a small correction determined by calibration. The range of the instrument without extrapolation is from 12% to 97% relative humidity. Over the range in relative humidity from 0% to 100% the voltage output ranges from 0 to 0.1 volt.

The calibration of the instrument is readily performed by immersion of the

probe into the vapor region above various saturated salt solutions [12]. The ones used here are lithium chloride (12%), magnesium chloride (33%), and potassium sulfate (97%). The calibration checks performed by some of the participating laboratories were not made by this method but by direct comparison to other similarly calibrated instruments. This form of calibration check holds true for the temperature and pressure instrumentation as well.

F. Multiplexer

The multimeter used to digitize the analog signals for acquisition by the computer has only a single input port and requires an external switching network for multiple inputs. Cutkosky designed a computer-controlled multiplexer. With his cooperation, we constructed a replica for this project. The multiplexer permits the computer software to control the presentation of the analog signals for acquisition by the computer. The multiplexer functioned well during the entire round-robin sequence. However, a malfunction would not necessarily be noticeable by an operator and a self-checking software routine verified that each of the four switching relays would trip and the contact resistance was acceptable before operation began in a participant laboratory.

In testing the multiplexer the operator had to become an active participant in the process, first loading into the computer and then running the appropriate program. The computer would cycle all of the relays, leave the switches in the open position and then, using the multimeter in the resistance mode, verify each switch is open. If all is well the computer instructs the operator to short each switch input with a jumper wire and then proceeds to see if the closed switch contact resistance is acceptable.

The cardinal virtues of the multiplexer are summarized below:

- o IEEE-488 control bus.
- o 128 relays are controllable.
- o No applied power after tripping relay.
- o Low thermal EMF.
- o High leakage resistance.
- o Zero offset current.
- o Gold switch contacts for low level signals.

G. Multimeter

The signal outputs from the instrumentation discussed thus far are DC voltages that range from a few millivolts to no more than 10 volts. To measure these voltages and to convert them into digital information useful to the computer we selected a programmable multimeter with an IEEE-488 interface. Thus we can control the meter setup, data collection, and transfer to the computer by software commands from the computer. As we have already discussed, software control of the multiplexer determines what signal is presented to the multimeter.

There are many features provided in the meter that we did not make use of however, the cardinal features we took advantage of are listed below:

- o Five and one-half digit resolution
- o Auto ranging
- o Remote programing
- o IEEE-488 interface
- o External trigger
- o High input impedance

H. Computer

The Hewlett-Packard 85 personal computer was chosen for this project for several reasons. Below is a list of properties that we felt were desirable:

- o Self-contained
- o Rugged
- o Highly portable
- o Magnetic tape cartridge mass storage
- o User-friendly interfacing
- o Small size

The HP-85 computer contains in one frame a 5-inch cathode ray tube display, keyboard, tape cassette, and a 32-column printer. To increase the computer's capability we added four read only memory units (ROM's), a 16k memory expansion module, an IEEE interface unit and an external printer. The added ROM's are:

- o External input-output
- o Matrix arithmetic enhancement
- o Printer/Plotter enhancement
- o Advanced programming

With these enhancements the computer is able to control the measurement, permit easy editing of data files, reduce the data, analyze the results, and output the desired information to an 80-column printer.

An essential part of the computer is the traveling case obtained from the manufacturer. We believe this case and the package shipping container make it possible to repeatedly ship the computer and the other equipment without damage. The shipping case is made of molded plastic and supports the computer on antishock mountings while in transit. The larger external shipping container is designed to withstand the abuse received from 100 air freight shipments and meets the Air Transportation Association Specifications 300 Category 1. Internally the computer and the other instruments are each surrounded by Dow #220 polyethylene foam sheets that provides excellent shock protection. The outside shape of the shipping container is cubic with a distance of approximately 1 meter separating the cube faces.

I. Computer Software

The three programs used in this package have been recorded onto two magnetic tape cartridges as shown in table 3. The first cartridge contains two programs which are used to test the integrity of the data acquisition system (i.e., the multiplexer, multimeter, thermometer, barometer, and relative humidity indicator). The second tape cartridge stores one main program, its subprograms, and data files which are used to record observations, analyze data, and print the results of a complete set of weighings. A sample printout of these results can be found in table 1 of appendix III. In addition, appendices I and II will describe the content and function of each tape cartridge while providing the reader with a printout of the programs found in table 3.

Results

The participants were given written instructions regarding the proper assembly, test and operation of the equipment. Certain aspects of the measurement preparation such as weight soaking time, balance warmup, and adjustment of the fixed heat source were dependent on the operator following the instructions carefully. Three laboratories (Sandia, Mound, and 3M) were asked to check the calibration of the thermistor thermometer, barometer and humidity indicator and inform us if they found any significant change. Nothing irregular was detected. However, we caused a large shift in the barometer calibration after completion of the measurements while performing a recalibration at the NIST. Weight H2 has a very large buoyancy correction and the calculated mass is sensitive to barometer error. Our data make us confident the shift occurred after use at the Mound laboratory and not beforehand. Recalibration of the other instruments was straightforward and no significant differences were found.

On return of the package to the NIST we repeated the measurements again. The data clearly indicate that weight H2 lost 1/2 mg of mass after calibration at Mound Laboratories and before final calibration at the NIST. This loss correlates with an obvious deep scratch on the weight surface. Although we show the data for completeness we do not use it in the analysis. Other discrepancies in the data occurred as well. Sandia Laboratories was unable to weigh kilogram TE because of its height and a previous modification of their balance pan. There are two sets of Hawaiian data, the second set is a recomputation of the first set after an NIST recalibration of their kilogram standards. Colorado violated a procedural rule that limits the difference between the standard and the unknown weight to no more than 1/4 the mass of the sensitivity weight. The violation occurred on the A1-S2 weight pair and we omit the data. For completeness we present the data of our first round-robin (RR#1) along with the second one (RR#2) here. Tables 4, 5, and 6 present the first and second round-robin data in tabular form. Clearly, we have eliminated the difficulty the state laboratories had with the buoyancy correction for weight H2 during round-robin #1. We have improved the standard deviations for each weight that was included earlier and note that weight TE (not included in RR#1) appears to behave as the others. The t-test was failed

Table 3. Catalog of programs stored on the preliminary and primary tape cartridges.

Preliminary tape cartridge:

NAME	TYPE	BYTES	RECS	FILE
TEST	PROG	256	15	1
HIST	PROG	256	7	2

Primary tape cartridge:

NAME	TYPE	BYTES	RECS	FILE
Autost	PROG	256	5	1
INPUT	PROG	256	12	2
MEASUR	PROG	256	13	3
EDITD	PROG	256	23	4
ANALY	PROG	256	26	5
SUBRHO	PROG	256	4	6
PRINTM	PROG	256	4	7
A MAIN	DATA	80	12	8
B MAIN	DATA	80	12	9
C MAIN	DATA	80	12	10
TRANSD	PROG	256	9	11

in every instance for the New York data during the second round-robin and indicates the state kilograms need recalibration. This information is a useful byproduct of the round-robin.

Hawaii was the second participating laboratory to receive the package from NIST. We noted that based on our experience we thought the Hawaiian kilograms appeared to need recalibration. This was accomplished by us and explains the recalculation of the Hawaiian data. We found the mass of both kilograms in the uncleaned condition, as they were received from Hawaii, to be virtually unchanged but in preparing for an after cleaning calibration we routinely measured their densities. Recalculation of the data with the measured densities brings the Hawaiian measurements in agreement with the others. Recalibration of the Hawaii kilograms after cleaning showed a mass loss of about 150 micrograms per kilogram (i.e., the dirt compensated for the wear). More discussion on the matter of artifact density occurs later in the text.

In performing the measurements reported here, each participating laboratory used their pair of calibrated standard kilogram to make mass assignments to the artifacts in the round-robin package. The assigned mass values of these standard kilogram pairs can be offset from the others by 100 micrograms or more based on their calibration uncertainty statements. These offsets then manifest themselves in the subsequent artifact calibration. We have removed these offsets, for the purpose of comparison, before performing additional analysis by normalizing the A1 mass values and then adjusting the other artifact data accordingly. Kilogram A1 was chosen to normalize the data because its characteristics are similar to the participants' standard kilograms. The adjustment factor for each participant was obtained by subtracting their averaged A1 values from the grand average of all the A1 values. This factor was then added algebraically to the assigned masses of the remaining kilograms. We refer to these data as the adjusted data and show it in table 7.

For RR#2 the adjustment factor never exceeds 100 micrograms and is a negative 100, 60, and 56 micrograms respectively for the states of Nevada, Utah, and Colorado -- all high-altitude locations. We shall speculate on the significance of this finding later in the text. Comparing the standard deviations of these tables to those given in tables 4, 5, and 6 we see an across-the-board improvement in the data for both round-robins.

We solicited the assistance of the Statistical Engineering Division to perform additional statistical tests on the data. In particular we wanted to look for correlation between the artifact masses and the air density parameters of barometric pressure, air temperature, and relative humidity. In addition we wanted to look for correlation with mass and the difference between the initial temperature and the final temperature found inside the balance during the weighing period. We refer to this difference as ΔT . Kilogram TE was specifically designed to be sensitive to this temperature difference [8]. These data are given in table 8 and the analysis in table 9. We quote the summary here:

Table 4. Assigned mass values for weight pair A1-S2, given for both round-robins.

Round-robin #1:		Kilogram Designation		Obs S.D. (μ g)	F Ratio ¹	Check Std	t- Test ²
Lab	Date	A1 1000 (g)	S2 999 (g)			Accept -Obs (μ g)	
NIST ³	4/82-4/83	.009132	.988553	-	-	-	-
New York	5/82	.009295	.988625	21	0.26	-86	3.04 ⁴
California	6/82	.009301	.988789	47	0.61	-4	0.09
Nevada	8/82	.009070	.988476	79	1.75	-71	-1.68
Sandia	11/82	.009102	.988739	61	2.74	14	-2.63
Colorado	2/83	.009068	.988602	54	1.72	4	-0.89
Utah	3/83	.009239	.988689	34	0.32	-38	0.90
\bar{X} =		.009172	.988782				
S.D. =		0.000103	0.000377				

Round-robin #2:		Kilogram Designation		Obs S.D. (μ g)	F Ratio	Accept Diff	t- Test
Lab	Date	A1 1000 (g)	S2 999 (g)			-Obs Diff (μ g)	
NIST #1	3/86	.009123	.988540	52	2.21	5	-0.21
NIST #2	3/86	.009191	.988539	25	0.53	-66	2.67
California	6/86	.009177	.988578	28	1.98	37	-2.61
Hawaii #1 ⁶	7/86	.009348 ⁵	.988990 ⁵	33	0.48	34	-1.00
Hawaii #2 ⁷	7/86	.009204	.988846	33	0.48	4	-0.12
3M Corp	8/86	.009173	.988625	82	3.59	37	-1.21
Nevada	9/86	.009270	.988804	25	0.12	-6	0.12
Utah	10/86	.009234	.988823	48	0.54	73	-1.59
Colorado ⁸	11/86	.009279 ⁵	.989044 ⁵	64	0.63	82	-1.43
Sandia	1/87	.009142	.988498	37	1.35	23	-1.01
New York	3/87	.009163	.988804	20	0.25	106	-3.74 ⁴
Mound	5/87	.009130	.988538	22	0.64	-4	0.20
NIST #3	8/87	.009126	.988390	44	1.56	-66	2.68
\bar{X} =		.009176	.988635				
S.D. =		0.000047	0.000157				

1. Critical value is 3.79 for probability = 0.01.

2. Critical value is 3.

3. NIST data is an average of seven measurements.

4. Failed t-test.

5. Data not included in the analysis.

6. Data reduction before recalibrating Hawaii standards.

7. Data reduction after recalibrating Hawaii standards.

8. Added weight rule violation.

Table 5. Assigned mass values for weight pair R2-H2, given for both round-robins.

Round-robin #1:

Lab	Date	Kilogram Designation		Obs S.D. (μ g)	F Ratio	Check Std	t- Test
		R2 999 (g)	H2 1000 (g)			Accept -Obs (μ g)	
NIST	4/82-4/83	.994503	.003880	-	-	-	-
New York	5/82	.994655	.005980	35	0.35	-62	1.45
California	6/82	.994742	.004544	42	0.49	-31	0.73
Nevada	8/82	.994471	.004643	65	1.16	84	-1.97
Sandia	11/82	.994436	.004143	60	2.85	2	-1.59
Colorado	2/83	.994580	.004642	65	0.82	-38	0.08
Utah	3/83	.994788	.003816	54	0.16	5	-0.12
\bar{X} =		.994596	.004521				
S.D. =		0.000136	0.000731				

Round-robin #2:

Lab	Date	Kilogram Designation		Obs S.D. (μ g)	F Ratio	Accept Diff	t- Test
		R2 999 (g)	H2 1000 (g)			-Obs Diff (μ g)	
NIST #1	3/86	.994417	.003937	33	0.88	-5	0.19
NIST #2	3/86	.994429	.003944	29	0.67	-68	2.77
California	6/86	.994453	.003973	24	1.43	16	-1.16
Hawaii #1	7/86	.994661 ¹	.004237 ¹	93	3.73	33	-0.98
Hawaii #2	7/86	.994517	.004093	93	3.73	3	-0.10
3M Corp	8/86	.994537	.004058	45	1.06	-30	0.99
Nevada	9/86	.994677	.004137	78	1.15	14	-0.28
Utah	10/86	.994708	.004094	46	0.50	99	-2.16
Colorado	11/86	.994804	.004384	71	0.76	117	-2.04
Sandia	1/87	.994476	.004037	56	3.08	-44	1.95
New York	3/87	.994566	.004109	34	0.70	102	-3.59 ²
Mound	5/87	.994494	.003949	49	3.32	21	-1.12
NIST #3	8/87	.994368	.003482 ^{1,3}	49	1.96	-57	2.30
NIST #4	8/87	.994390 ¹	.003504 ^{1,3}	62	3.18	-126	5.07 ²
\bar{X} =		.994537	.004065				
S.D. =		0.000131	0.000128				

1. Data not used in the analysis.
2. Failed t-test.
3. Damaged weight (see text).

Table 6. Assigned mass values for weight pair Al-TE. The pair Al-TE was not used in round-robin #1.

Round-robin #2:

Lab	Date	Kilogram Designation		Obs S.D. (μ g)	F Ratio	Accept Diff	t- Test
		Al 1000 (g)	TE 1000 (g)			-Obs Diff (μ g)	
NIST #1	3/86	.009195	.002105	59	2.81	-34	1.36
NIST #2	3/86	.009096	.002082	53	2.32	-38	1.54
California	6/86	.009226	.002277	46	5.28	-33	2.37
Hawaii #1	7/86	.009261 ¹	.002520 ¹	55	1.33	101	-2.98
Hawaii #2	7/86	.009117	.002376	55	1.33	71	-2.10
3M Corp	8/86	.009211	.002373	59	1.83	-5	0.18
Nevada	9/86	.009285	.002457	70	0.92	-34	0.65
Utah	10/86	.009242	.002389	74	1.29	-32	0.69
Colorado	11/86	.009188	.002519	23	0.84	69	-1.20
Sandia ¹	1/87	-	-	-	-	-	-
New York	3/87	.009187	.002273	34	0.71	86	3.05 ²
Mound	5/87	.009127	.002301	14	0.28	20	-1.05
NIST #3	8/87	.009114	.002059	57	2.65	-55	2.23
\bar{X} =		.009181	.002292				
S.D. =		0.000060	0.000154				

-
1. Modified balance pan cannot accommodate extraordinary height of weight TE.
 2. Failed t-test.

Table 7. Assigned mass values of weights S2, R2, H2, and TE normalized to weight A1.

Round-robin #1:

		Kilogram Designation			Adjust Factor A1 - A1 (μ g)
Lab	Date	S2 999 (g)	R2 999 (g)	H2 1000 (g)	
NIST	4/82-4/83	.988593	.994543	.003920	40
New York	5/82	.988502	.994532	.005857	-123
California	6/82	.988660	.994613	.004415	-129
Nevada	8/82	.988578	.994573	.004745	102
Sandia	11/82	.988809	.994506	.004213	70
Colorado	2/83	.988706	.994684	.004746	104
Utah	3/83	.988622	.994721	.003749	-67
\bar{X} =		.988639	.994596	.004521	
S.D. =		0.000099	0.000081	0.000701	

Round-robin #2:

		Kilogram Designation				Adjust Factor A1 - A1 (μ g)
Lab	Date	S2 999 (g)	R2 999 (g)	H2 1000 (g)	TE 1000 (g)	
NIST #1	3/86	.988559	.994436	.003956	.002124	19
NIST #2	3/86	.988573	.994462	.003978	.002116	34
California	6/86	.988554	.994429	.003949	.002253	-24
Hawaii	7/86	.988864	.994535	.004111	.002394	18
3M Corp	8/86	.988611	.994523	.004044	.002359	-14
Nevada	9/86	.988704	.994577	.004037	.002357	-100
Utah	10/86	.988763	.994648	.004034	.002329	-60
Colorado	11/86	1	.994748	.004328	.002463	-56
Sandia	1/87	.988534	.994512	.004073	2	36
New York	3/87	.988807	.994569	.004112	.002276	3
Mound	5/87	.988588	.994544	.003999	.002351	50
NIST #3	8/87	.988448	.994426	3	.002117	58
\bar{X} =		.988637	.994534	.004056	.002285	
S.D. =		0.000129	0.000095	0.000106	0.000120	

1. Added weight violation.
2. Balance modification.
3. Damaged weight.

Table 8. Correlation test data.

Kilogram	Adj. Mass (g)	Temperature (° C)	Pressure (mm Hg)	Humidity (%)	Δ Temp (° C)
S2	999.988559	24.08	761.25	32.7	-0.01
	999.988573	24.24	752.44	27.8	0.02
	999.988554	23.30	759.58	38.4	0.00
	999.988864	23.36	762.30	23.8	0.01
	999.988611	22.67	737.26	33.5	-0.01
	999.988704	24.40	646.41	18.8	0.09
	999.988763	23.74	656.26	26.8	0.10
	999.988534	20.87	627.17	38.1	-0.02
	999.988807	21.20	762.44	40.6	0.18
	999.988588	22.96	743.61	33.4	0.04
	999.988448	24.20	748.63	53.5	0.00
R2	999.994436	24.27	756.77	34.0	0.02
	999.994462	24.05	759.64	27.5	0.03
	999.994429	23.27	759.13	38.4	0.00
	999.994535	23.39	763.10	24.0	-0.06
	999.994523	22.59	743.66	32.8	0.00
	999.994577	24.14	645.08	17.0	0.15
	999.994648	23.12	655.87	32.7	0.07
	999.994748	22.37	619.10	32.3	-0.02
	999.994512	20.84	626.63	37.3	-0.10
	999.994569	22.79	751.60	33.6	0.01
	999.994544	22.97	743.44	33.5	0.02
	999.994426	24.24	754.14	49.2	0.01
H2	1000.003956	24.27	756.77	34.0	0.02
	1000.003978	24.05	759.64	27.5	0.03
	1000.003949	23.27	759.13	38.4	0.00
	1000.004111	23.39	763.10	24.0	-0.06
	1000.004044	22.59	743.66	32.8	0.00
	1000.004037	24.14	645.08	17.0	0.15
	1000.004034	23.12	655.87	32.7	0.07
	1000.004328	22.37	619.10	32.3	-0.02
	1000.004073	20.84	626.63	37.3	-0.10
	1000.004112	22.79	751.60	33.6	0.01
	1000.003999	22.97	743.44	33.5	0.02
TE	1000.002124	24.29	741.27	41.4	-0.01
	1000.002116	24.21	753.34	30.7	0.02
	1000.002253	23.25	759.50	35.6	0.06
	1000.002394	23.52	762.45	23.3	0.13
	1000.002359	22.61	745.70	33.2	-0.01
	1000.002357	25.47	641.38	16.1	0.15
	1000.002329	23.58	657.15	28.1	0.01
	1000.002463	21.68	630.45	30.6	0.13
	1000.002276	22.96	744.21	36.6	0.09
	1000.002351	23.01	743.24	33.1	0.03
	1000.002117	24.34	755.54	49.7	0.00

Table 9. Correlation of Mass with Environmental Parameters.

Kilogram Designation	Quantity	Correlation Coefficient	Correlation Coefficient 95% Confidence Interval		
S2	Temperature	-0.14	-0.68	+0.50	
	Pressure	-0.03	-0.62	+0.58	
	Humidity	-0.58	-0.88	+0.03	
	Δ Temperature	+0.64	+0.06	+0.90	*
R2	Temperature	-0.38	-0.78	+0.25	
	Pressure	-0.70	-0.91	-0.21	*
	Humidity	-0.36	-0.77	+0.27	
	Δ Temperature	+0.12	-0.49	+0.65	
H2	Temperature	-0.44	-0.82	+0.22	
	Pressure	-0.54	-0.86	+0.08	
	Humidity	-0.08	-0.65	+0.54	
	Δ Temperature	-0.29	-0.76	+0.38	
TE	Temperature	-0.53	-0.86	+0.10	
	Pressure	-0.53	-0.86	+0.10	
	Humidity	-0.64	-0.90	-0.06	*
	Δ Temperature	+0.62	+0.03	+0.89	*

Nevada, Utah, and Colorado excluded:

S2	Temperature	-0.31	-0.81	+0.44	
	Pressure	+0.04	-0.64	+0.68	
	Humidity	-0.53	-0.88	+0.21	
	Δ Temperature	+0.57	-0.15	+0.89	
R2	Temperature	-0.55	-0.89	+0.18	
	Pressure	-0.19	-0.76	+0.54	
	Humidity	-0.50	-0.87	+0.25	
	Δ Temperature	-0.28	-0.80	+0.47	
H2	Temperature	-0.51	-0.89	+0.30	
	Pressure	-0.28	-0.82	+0.53	
	Humidity	-0.34	-0.84	+0.48	
	Δ Temperature	-0.57	-0.91	+0.22	
TE	Temperature	-0.84	-0.97	-0.33	*
	Pressure	+0.07	-0.67	+0.74	
	Humidity	-0.67	-0.93	+0.07	
	Δ Temperature	+0.53	-0.27	+0.90	

* Indicates significant correlation (95% confidence interval does not overlap zero).

Table 10. - Measurement series temperature change.

T_i = initial temperature

T_f = final temperature

$\Delta T = T_f - T_i$

Round-robin #2:

Lab	Temperature Difference (ΔT) (° C)		
	A1-S2	R2-H2	A1-TE
NIST #1	-0.01	0.02	-0.01
NIST #2	0.02	0.03	0.02
CA	0.00	0.00	0.06
Hawaii	0.01	-0.06	0.13
3M Corp	-0.01	0.00	-0.01
Nevada	0.09	0.15	0.15
Utah	0.10	0.07	0.01
Colorado	-	-0.02	0.13
Sandia	-0.02	-0.10	-
New York	0.18	0.01	0.09
Mound	0.04	0.02	0.03
NIST #3	0.00	0.01	0.00
NIST #4	-	0.05	-

"Note that on an individual basis the correlation is significant at the 0.05 level in only four cases, as indicated by a "*." However, when considering the data for the four kilograms as whole, there are strong indications that

- 1) mass is negatively correlated with temperature, pressure, and humidity,
- 2) mass is positively correlated with delta temperature."

Only kilogram R2 has significant correlation with barometric pressure and therefore we suspected the weight volume was in error. However, a new volume determination agreed with the former one to within 1 part in 10^4 , a null result. We expected correlation with kilogram TE due to thermal effects [8] but not the general correlation of all the kilograms with all of the parameters, particularly relative humidity. Our data yields a negative correlation where a positive correlation is known [13]. We recalled during the volume measurement on Hawaii's kilogram that the assigned volumes to state kilograms were assumed and not measured at the time of their calibration. This occurred during the 1970's and before, when volume measurements were not economical. We performed a second analysis but omitted the data from Colorado, Nevada, and Utah. This analysis is summarized in the lower part of table 9 and indicates only one significant correlation, kilogram TE with delta T. The kilograms and the observed delta T are given in table 10. There remains a suggestion of a negative correlation of mass with temperature and relative humidity that we cannot explain other than possible instrumentation errors. However, the resulting mass errors are well below the noise level of these measurements or any practical routine mass measurement. Fig. 4 is a typical graphic illustration of both sets of analysis, we only show the data for kilogram R2 where earlier we had considered a volume error. In Fig. 5 we have denoted the weight pairing of kilogram A1 at the time of weighing. It is gratifying to note that the mass assigned to kilogram A1 is not correlated with the unusual features of kilograms S2 and TE. Other aspects of Fig. 5 are discussed later in the text.

Conclusions


First, let us consider only the improvement we observe in our ability to weigh kilogram H2. The volume difference between H2 and the participants standards' is about two orders of magnitude larger than the other kilogram combinations. We note that in round-robin #1 the weighing data was unacceptable while that of round-robin #2 is as good as can be expected [6] and not significantly different from the behavior of the other kilograms. This certainly indicates that the improvements in instrumentation and thermal conditioning are working very well.

Secondly, the mass of kilogram TE is the only mass that has a significant correlation with delta T. The height of TE is about three times that of a standard kilogram, similar to A1, whose height is approximately equal to its diameter. Furthermore we note that the surface area of kilogram S2 is about 4.5 times larger than a standard kilogram and no significant correlation with any parameter was detected. These data suggest that we

Figure 4. Scatter plots of correlated and uncorrelated data for kilogram R2.

**Correlation Plots of Mass Values Assigned to Kilogram R2
Versus Temperature, Pressure, Relative Humidity, and ΔT**

With Nevada, Utah, and Colorado

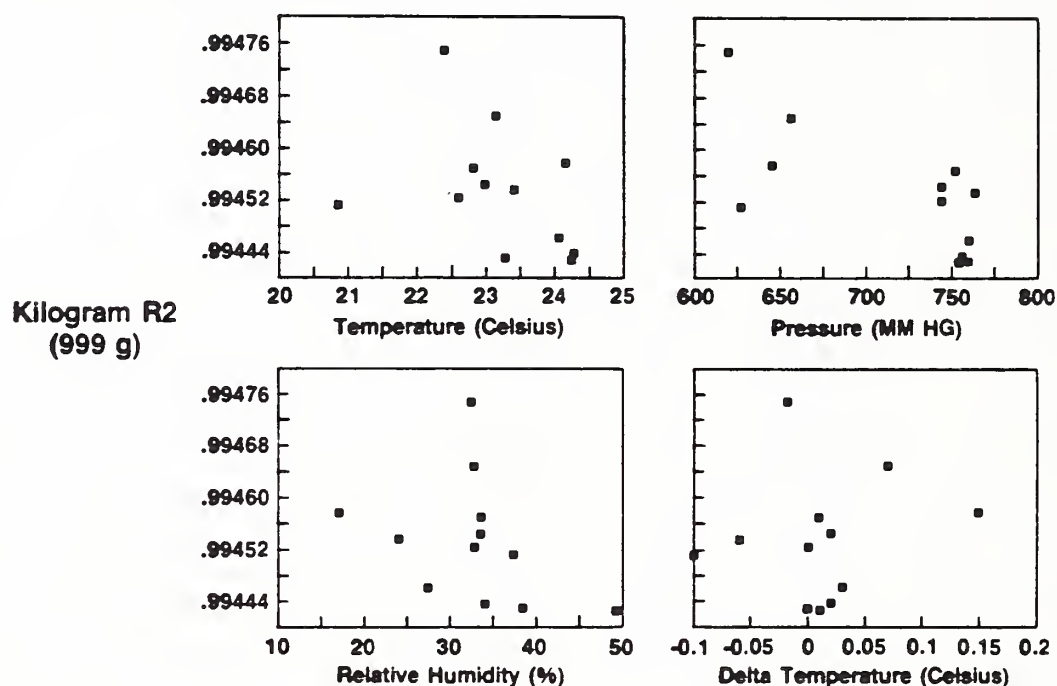

Without Nevada, Utah, and Colorado


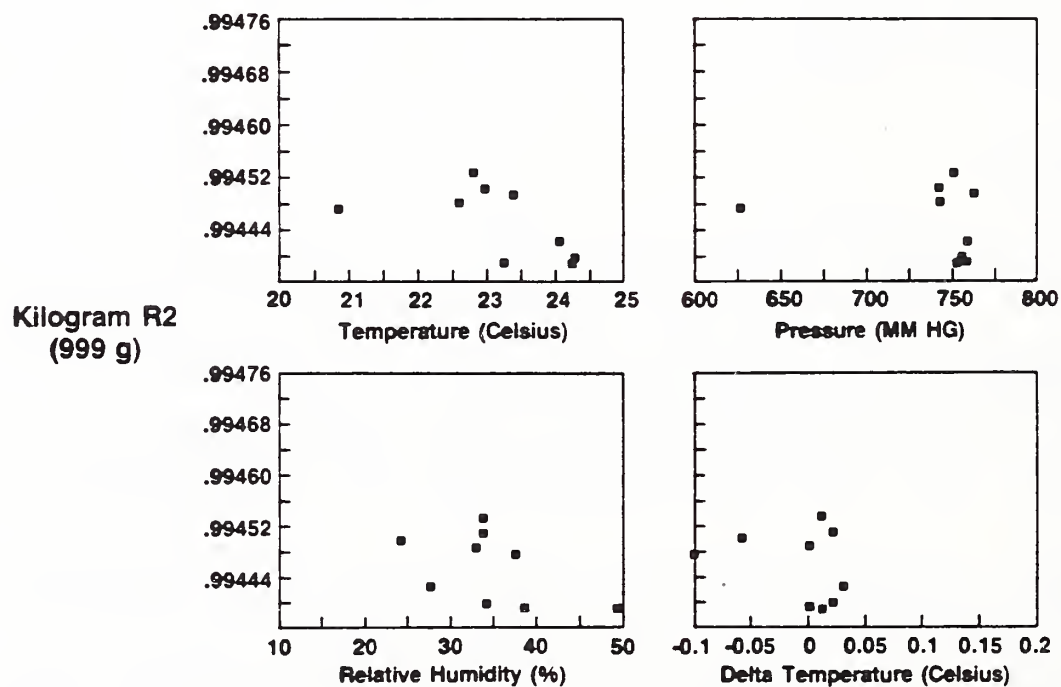


Figure 5. These data illustrate the effects of assumed density on the assigned mass of kilogram A1.

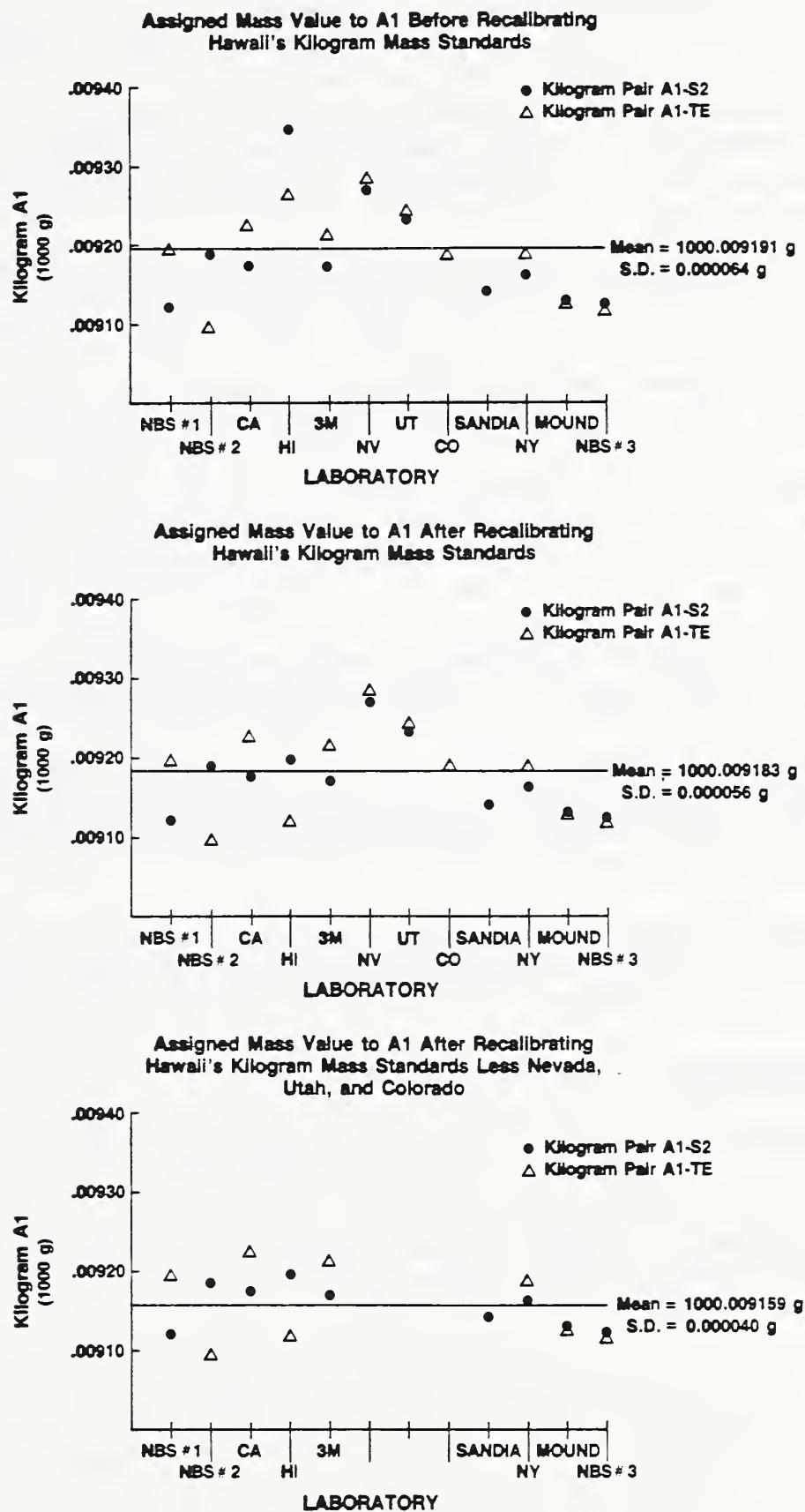

could extend our measurement process to include a summation of 1 kilogram weights (500g, 300g, 200g) that have less exaggerated characteristics, without significant systematic errors. This is an important consideration in the calibration of a set of weights if the uncertainty estimates are to be meaningful.

We also note the behavior of kilogram A1 with the data plots of Fig. 5. Although the kilogram pairing is denoted by dots and triangles it has no significance in this discussion. The above plot includes the Hawaiian point based on a faulty volume assumption and the middle plot with the correct volume. In the lower plot we have removed the data from the states of Colorado, Nevada, and Utah, each with assumed volumes for their standards and used in laboratories at higher elevations and therefore lower air densities than the other laboratories. The standard deviation is not only improved but is very good for routine mass calibrations and we speculate the middle plot would likewise be improved with measured weight volumes used in the data reduction. The New York kilograms, like nearly all state kilograms, have not had their volumes measured and it is probably fortuitous that our data do not reflect a faulty volume assumption. That is, we know the actual material density is somewhat lower than the assumed 8.0 g/cm^3 and that the weights were calibrated and subsequently used at about the same air density thus masking the error. Furthermore the statistical check indicates the New York Kilograms are in need of recalibration and may have some compensating effect for the error. We believe we can use an improved version of this package to calibrate mass standards in the recipient laboratory provided their kilograms have had their densities adequately determined beforehand. Any future effort would be directed toward a complete weight set calibration. Communication with the participants revealed that most, if not all, were pleased with the concept of the package and with the learning experience arising from performing the measurements. Indeed, we believe that NIST successfully transferred near state-of-the-art mass measurement technology to the individuals that participated in the program.

After several years of experience interacting with the state laboratories we believe a program that sponsors the acquisition of computers, software, and electronic instrumentation is needed along with the required support.

Acknowledgements

We appreciate the funding support provided over the past few years by both Albert Tholen, Chief of the Office of Weights and Measures, and Joe Simmons, our former Division Chief. We are also grateful to Charles Reeve of the Statistical Engineering Division who always provides his statistical advice with good humor.

REFERENCES:

1. H. A. Bowman, W. Gallagher, and R. M. Schoonover, "The development of a working density standard," Instrument Society of America, 20th annual ISA conference and exhibit, Oct. 4-7, 1965.
2. P. E. Pontius, "Mass Measurements: A study of anomalies," Science, vol. 190, pp. 379-380, 1975.
3. R. M. Schoonover and J. E. Taylor, "Some recent developments at NBS in mass measurements," IEEE Transactions on Instrumentation and Measurement, vol. IM-35, no. 4, part 1, Dec. 1986.
4. H. A. Bowman and R. M. Schoonover, "Procedure for high precision density determination by hydrostatic weighing," J. Res., vol. 71C, pp. 179-198, July - Aug. 1967.
5. F. E. Jones, "The air density equation and the transfer of the mass unit," J. Res., vol. 10, pp. 419-427, Sept. - Oct. 1978.
6. R. M. Schoonover, R. S. Davis, R. G. Driver, and V. E. Bower, "A practical test of the air density equation in standards laboratories at differing altitudes," J. Res., vol. 85, no. 1, pp 27-38, Jan. - Feb. 1980.
7. R. M. Schoonover, "Advanced mass calibration in state laboratories," Rep. NBSIR 83-2752, pp. 1-9, Aug. 1983.
8. R. M. Schoonover and J. Keller, "A surface dependent thermal effect in mass calibration," in "Report of the 68th National Conference on Weights and Measures 1983," NBS Special Pub. 663, pp. 39-57, 1983.
9. M. E. Cage and R. S. Davis, "An analysis of read-out perturbations seen on an analytical balance with a swinging pan," J. Res. vol. 87, no. 1, Jan. - Feb. 1982.
10. H. A. Bowman and H. E. Almer, "Minimization of the arrestment error in one-pan, two-knife balance systems," J. Res. Nat. Bur. Stand. (U.S.), 67C (Eng. and Instr.), No. 3, 227-235 (July-Sept. 1963).
11. R. D. Cutkosky and R. S. Davis, "Simple control circuit for temperature regulation and other bridge applications," Rev. Sci. Instrum., vol. 52, no. 9, pp. 1403-1405, Sept. 1981.
12. A. Wexler and S. Hasegawa, "Relative humidity-temperature relationships of some saturated salt solutions in the temperature range 0° to 50° C," J. Res., vol. 53, no. 1, July 1954.

13. M. Kochsiek, "Measurement of Water Adsorption on Metal Surfaces," *Metrologia* 18, 153-159, 1982.

Appendices

Appendix I.

Preliminary tape cartridge.....32

Appendix II.

Primary tape cartridge.....38

Appendix III.

Printout of results.....57

Appendix I.

In conjunction with the computer, the preliminary tape cartridge is used as a debugging tool during the set up of the data acquisition system. This auxiliary tape, contains the programs TEST and HIST.

TEST is a program that functions as a continuity tester and is used to check the opening and closing of the multiplexer relays before connecting the multiplexer to the remaining components of the data acquisition system.

The program HIST will allow the user to automatically measure temperature, pressure, and relative humidity if the data acquisition system is functioning properly. In addition to using HIST as a diagnostic tool, each laboratory can periodically monitor room temperature and choose the most quiescent period of the day in which to make measurements.

```

1000 REM ***** TEST *****
1010 REM
1020 REM
1030 REM THIS PROGRAM IS A TEST
1040 REM FOR CONTINUITY OF THE
1050 REM 3 RELAYS USED IN THE
1060 REM SECOND ROUND ROBIN.
1070 REM
1080 REM
1090 REM *** INITIALIZE ***
1100 REM
1110 PRINTER IS 2
1120 OPTION BASE 1
1130 INTEGER P2(7,5),P3(7,5),P4(7,5),M(7,5),T(7,5)
1140 INTEGER I,J,N,K
1150 REAL R(7,5),V,O,C
1160 O=9000000000
1170 C=.4
1180 N=1
1190 CLEAR
1200 DISP USING "2/" ;
1210 DISP "1) CHECK ALL CONNECTIONS BEFORE"
1220 DISP "      RUNNING THIS TEST."
1230 DISP USING "1/" ;
1240 DISP "2) THE PROGRAM IS APPROX.          5 MINUTES LONG."
1250 DISP USING "1/" ;
1260 DISP "3) RESULTS WILL BE DISPLAYED ON    THE SCREEN AND/OR OUTPUT"
1270 DISP "      TO THE THERMAL PRINTER."
1280 DISP USING "1/" ;
1290 DISP "4) PRESS THE CONT KEY TO RUN      THE TEST."
1300 PAUSE
1310 CLEAR
1320 REM
1330 REM *** P2,P3,P4 ***
1340 REM
1350 FOR J=1 TO N
1360 FOR I=1 TO 7
1370 P2(I,J)=0
1380 IF I=2 THEN P2(I,J)=1
1390 P3(I,J)=0
1400 IF I=4 THEN P3(I,J)=1
1410 P4(I,J)=0
1420 IF I=6 THEN P4(I,J)=1
1430 NEXT I
1440 NEXT J
1450 REM
1460 REM *** SET RELAYS ***
1470 REM
1480 FOR J=1 TO N
1490 FOR I=1 TO 7
1500 IF P2(I,J) OR P3(I,J) OR P4(I,J)=1 THEN 1540
1510 OUTPUT 710 ;"C1-1-2-3-4-5-6-7-8"
1520 GOSUB 1790
1530 GOTO 1730
1540 IF P2(I,J)=0 THEN 1580
1550 OUTPUT 710 ;"C1+2-3-4"
1560 GOSUB 1790
1570 GOTO 1730
1580 IF P3(I,J)=0 THEN 1620
1590 OUTPUT 710 ;"C1-2+3-4"

```

```

1600 GOSUB 1790
1610 GOTO 1730
1620 IF P4(I,J)=0 THEN 1660
1630 OUTPUT 710 ; "C1-2-3+4"
1640 GOSUB 1790
1650 GOTO 1730
1660 CLEAR
1670 DISP USING "4/" ;
1680 DISP "*****"
1690 DISP "PROGRAM ERROR!! ALL PORTS ARE OPEN."
1700 DISP "*****"
1710 BEEP @ BEEP @ BEEP @ BEEP @ BEEP
1720 PAUSE
1730 NEXT I
1740 NEXT J
1750 GOTO 1930
1760 REM
1770 REM *** READ METER ***
1780 REM
1790 WAIT 3000
1800 OUTPUT 724 ; "F3RAN4T2"
1810 R(I,J)=0
1820 FOR K=1 TO 10
1830 WAIT 300
1840 TRIGGER 724
1850 ENTER 724 ; V
1860 R(I,J)=R(I,J)+V
1870 NEXT K
1880 R(I,J)=R(I,J)/10
1890 RETURN
1900 REM
1910 REM *** TEST READING ***
1920 REM
1930 FOR J=1 TO N
1940 FOR I=1 TO 7
1950 IF R(I,J)>C THEN M(I,J)=0
1960 IF R(I,J)<O THEN M(I,J)=1
1970 IF R(I,J)>=C AND R(I,J)<=O THEN M(I,J)=10
1980 ON I GOTO 1990,2010,1990,2010,1990,2010,1990
1990 IF M(I,J)=0 THEN T(I,J)=0 ELSE T(I,J)=1
2000 GOTO 2020
2010 IF M(I,J)=1 THEN T(I,J)=0 ELSE T(I,J)=1
2020 IF R(I,J)>=C AND R(I,J)<=O THEN T(I,J)=1
2030 NEXT I
2040 NEXT J
2050 REM
2060 REM *** RESULTS ***
2070 REM
2080 FOR J=1 TO N
2090 FOR I=1 TO 7
2100 IF T(I,J)=1 THEN GOTO 2190
2110 NEXT I
2120 NEXT J
2130 CLEAR
2140 DISP USING "6/" ;
2150 DISP "THE MULTIPLEXER RELAYS ARE IN"
2160 DISP " GOOD WORKING ORDER."
2170 BEEP @ BEEP
2180 GOTO 2420
2190 CLEAR

```

```

2200 DISP USING "4/" ;
2210 DISP "WARNING!!!!!!!!!!!!!!"
2220 DISP USING "1/" ;
2230 DISP "THE MULTIPLEXER RELAYS HAVE          FAILED THE CONTINUITY TEST."
2240 DISP "PLEASE EXAMINE ALL ELECTRICAL    CONNECTIONS OR DETERMINE WHICH  RELAY
(S) IS FAULTY."
2250 BEEP @ BEEP @ BEEP @ BEEP @ BEEP
2260 PRINT USING "2/" ;
2270 FOR J=1 TO N .
2280 PRINT "N=";J
2290 FOR I=1 TO 7
2300 PRINT "      PORT   :   4 3 2   METER"
2310 PRINT USING 2320 ; P4(I,J),P3(I,J),P2(I,J),M(I,J)
2320 IMAGE 4X,"VALUE :",2X,3(D,X),3X,DD
2330 IF T(I,J)=1 THEN GOSUB 2390
2340 PRINT USING "2/" ;
2350 NEXT I
2360 PRINT USING "3/" ;
2370 NEXT J
2380 GOTO 2420
2390 PRINT USING 2400 ; R(I,J)
2400 IMAGE X,"V= ",K,X,"OHMS"
2410 RETURN
2420 END

```

```

1000 REM ***** HIST *****
1010 REM
1020 REM
1030 REM PROGRAM WILL AUTOREAD
1040 REM  ENVIROMENTAL DATA,
1050 REM  TIME, AND DATE.
1060 REM
1070 REM VALUES ARE OUTPUT TO
1080 REM  THE THERMAL PRINTER.
1090 REM
1100 REM TIME INTERVALS ARE
1110 REM  INPUT IN HOURS.
1120 REM
1130 REM
1140 PRINTER IS 2
1150 OPTION BASE 1
1160 REM *** SET DATE/TIME ***
1170 CLEAR
1180 DISP USING "3/" ;
1190 DISP "INPUT THE CURRENT DATE USING THE FOLLOWING FORMAT:"
1200 DISP USING "1/" ;
1210 DISP "          XX/XX/XXXX"
1220 DISP USING "1/" ;
1230 DISP " EXAMPLE: 04/05/1986"
1240 INPUT A$
1250 A=MDY(A$)-MDY("01/01/1985")
1260 CLEAR
1270 DISP USING "3/" ;
1280 DISP "INPUT THE CURRENT MILITARY   TIME USING THE FOLLOWING FORMAT:"
1290 DISP USING "1/" ;
1300 DISP "          HH:MM:SS"
1310 DISP USING "1/" ;
1320 DISP " EXAMPLE: 14:07:00"
1330 INPUT B$
1340 SETTIME HMS(B$),A
1350 REM *** SET TIMER ***
1360 CLEAR
1370 DISP USING "3/" ;
1380 DISP "INPUT TIME INTERVAL (HRS). "
1390 INPUT H1
1400 H1=3600000*H1
1410 CLEAR
1420 GOTO 1460
1430 ON TIMER# 1,H1 GOTO 1460
1440 GOTO 1440
1450 REM *** READ AMBIENT ***
1460 OUTPUT 724 ;"F1RAN5T2"
1470 OUTPUT 710 ;"C1-1+2-3-4-5-6-7-8"
1480 GOSUB 1600
1490 ENTER 724 ; T
1500 T=84.02531*T-.101
1510 OUTPUT 710 ;"C1-1-2+3-4-5-6-7-8"
1520 GOSUB 1600
1530 ENTER 724 ; P
1540 P=(30*P+800)/1.333224
1550 OUTPUT 710 ;"C1-1-2-3+4-5-6-7-8"
1560 GOSUB 1600
1570 ENTER 724 ; H
1580 H=H/.02-2.5
1590 GOTO 1640

```

```
1600 WAIT 2000
1610 TRIGGER 724
1620 RETURN
1630 REM *** PRINT DATA ***
1640 PRINT USING "2/" ;
1650 PRINT MDY$(DATE+MDY("01/01/1985"))
1660 PRINT TIMES$
1670 PRINT USING 1680 ; T
1680 IMAGE " T =",M3DZ.2D
1690 PRINT USING 1700 ; P
1700 IMAGE " P =",M3DZ.2D
1710 PRINT USING 1720 ; H
1720 IMAGE " H =",M3DZ.1D
1730 OFF TIMER# 1
1740 GOTO 1430
1750 END
```

Appendix II.

Each participating laboratory receives a duplicate copy of the primary tape cartridge with which weighings can be recorded, analyzed, and then returned to the NIST for a group comparison. This cartridge has one central program, seven subprograms, and three data files, table 3 in the main text. The three data files (A MAIN, B MAIN, and C MAIN) are storage areas in which all of the information associated with three sets of weighings (A1-S2, R2-H2, and A1-TE) are saved and modified by the various subprograms.

The main program "Autost" and its subprograms follow the logic outlined in Fig. 1. Whereby the user first selects a storage area (data file) followed by one of three processes (subprograms):

1. Input data using subprograms INPUT or MEASUR,
2. Revise the data file using subprogram EDITD, or
3. Compute and assign true mass values to the four kilograms associated with this data file using subprograms ANALY and SUBRHO.

The following is a brief description of the program files found on this tape cartridge:

Autost - The central program from which access to subprograms are obtained whenever necessary. As its name implies, the computer will automatically start executing commands from this program as soon as the user turns on power to the computer, provided the tape cartridge is in place.

INPUT - This subprogram allows the user to enter constants (the known standard deviation of the balance, volumetric coefficient of thermal expansion of the standard, etc.) from the keyboard into a data file.

MEASUR - Used to store, in a data file, balance readings and environmental conditions obtained in situ.

EDITD - An editor which allows the user to modify the data file whenever necessary.

ANALY - Contains the algorithm used to assign mass values to each kilogram and provide the user with a hard copy of results using the 80-column printer.

SUBRHO - Calculates an air density using temperature, pressure, and relative humidity values imported from subprogram ANALY.

PRINTM - A subprogram which will print the contents of a data file.

TRANSD - Although this program is not used during the weighings, it was written to be used as a mechanism through which data files can be transferred to and from tape cartridges.

Figure 1. Software flowchart.

1) LOAD AND RUN PROGRAM "Autost."

2) SELECT ONE OF THE THREE DATA FILES CORRESPONDING TO WEIGHING SERIES A1-S2, R2-H2, OR A1-TE.


3) USING SUBPROGRAM "INPUT," RECORD CONSTANT PARAMETERS OF THE WEIGHING SERIES INTO THE DATA FILE.

4) USING "MEASUR," RECORD BALANCE OBSERVATIONS AND AIR DENSITY READINGS INTO THE DATA FILE.

5) IF NECESSARY, MODIFY THE DATA FILE WITH "EDITD."

6) COMPUTE THE BOUYANCY CORRECTIONS AND MASS FOR EACH KILOGRAM USING "SUBRHO" AND "ANALY."

RESULTS ARE DISPLAYED WITH THE 80-COLUMN PRINTER AND "PRINTM" WILL LIST THE CONTENTS OF THE DATA FILE.

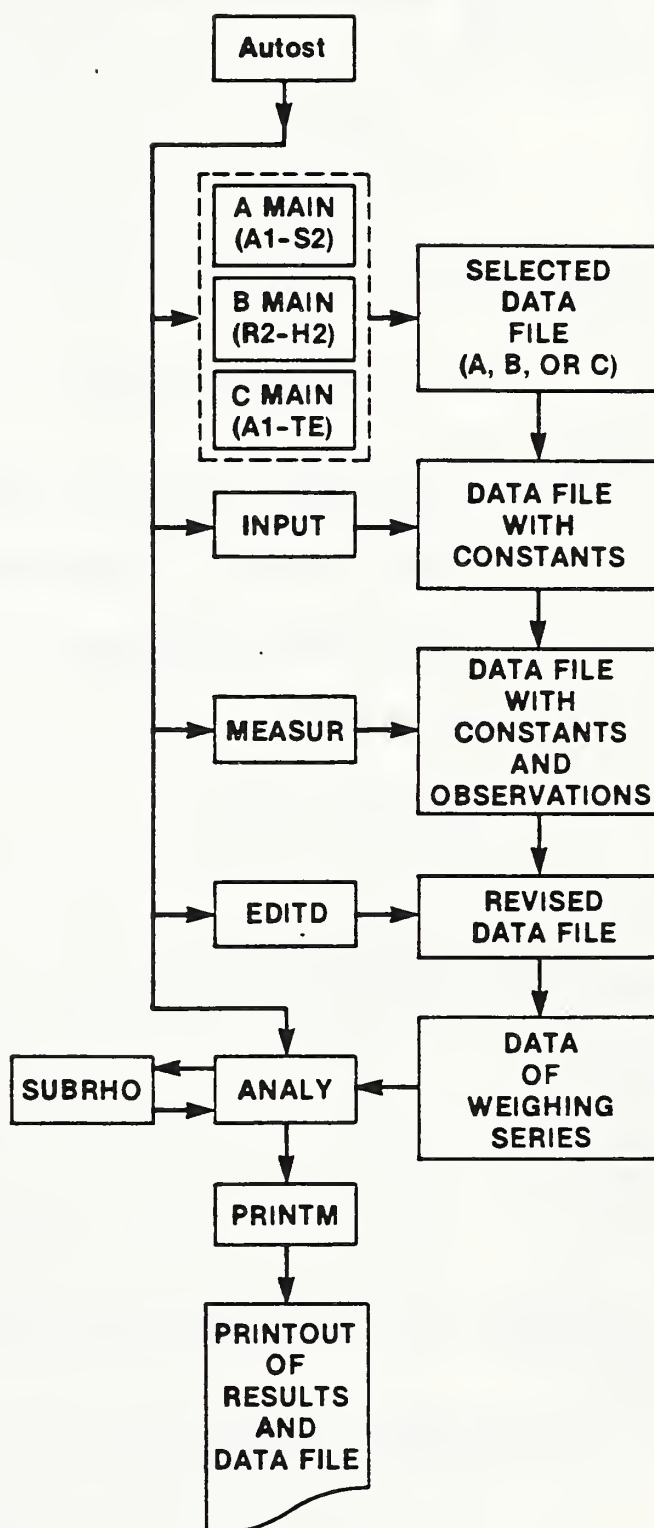

```

100 REM *MAIN PROGRAM
110 REM
120 REM *CONTROLLER FOR TRUE
130 REM MASS WEIGHINGS IN FOUR
140 REM /ONE'S SERIES.
150 REM
160 REM *WRITTEN BY J TAYLOR,
170 REM FOR OWM 3/85.
180 REM
190 CLEAR
200 PRINTER IS 701,80
210 DISP USING "2/"
220 DISP "TYPE IN MESUREMENT PAIR AS"
230 DISP "FOLLOWS :"
240 DISP USING "2/"
250 DISP "      1 = A1 S2"
260 DISP
270 DISP "      2 = R2 H2"
280 DISP
290 DISP "      3 = A1 TE"
300 DISP USING "2/"
310 INPUT X7
320 IF X7<1 OR X7>3 THEN GOTO 190
330 CLEAR
340 DISP USING "2/"
350 DISP "TYPE IN DESIRED PROGRAM FUNCTION AS FOLLOWS:"
360 DISP "
370 DISP "  1 = INPUT SYSTEM PARAMETERS."
380 DISP "
390 DISP "  2 = INPUT MEASUREMENTS."
400 DISP "
410 DISP "  3 = EDIT DATA FILES."
420 DISP "
430 DISP "  4 = ANALYZE DATA."
440 DISP
450 DISP "  5 = END"
460 INPUT X
470 IF X<1 OR X>5 THEN GOTO 330
480 ON X GOTO 490,510,530,550,590
490 CALL "INPUT" ( X7 )
500 GOTO 330
510 CALL "MEASUR" ( X7 )
520 GOTO 330
530 CALL "EDITD" ( X7 )
540 GOTO 330
550 CLEAR
560 CALL "ANALY" ( X7 )
570 CALL "PRINTM" ( X7 )
580 GOTO 330
590 CLEAR @ DISP USING "6/"
600 DISP "      RUN TERMINATED."
610 END

```

```

10 SUB "INPUT" (X7)
20 REM
30 REM *INPUTS DATA FROM KEY-
40 REM BOARD TO MAIN
50 REM DATA FILES.
60 REM
70 REM *WRITTEN BY J TAYLOR
80 REM FOR OWM 3/85.
90 REM
100 OPTION BASE 1
110 DIM H2$(20),H3$(8)
120 ON X7 GOTO 130,150,170
130 ASSIGN# 1 TO "A MAIN"
140 GOTO 180
150 ASSIGN# 1 TO "B MAIN"
160 GOTO 180
170 ASSIGN# 1 TO "C MAIN"
180 CLEAR
190 DISP
200 DISP "TYPE IN DESIRED PROGRAM FUNCTION AS FOLLOWS:"
210 DISP
220 DISP " 1 = INPUT HEADING/DATE."
230 DISP
240 DISP " 2 = INPUT CONSTANT DATA."
250 DISP
260 DISP " 3 = EXIT"
270 DISP
280 INPUT K
290 IF K<1 OR K>3 THEN GOTO 180
300 ON K GOTO 320,430,1270
310 REM *****HEADING/DATE*****
320 CLEAR
330 N=2
340 H2$=" "
350 DISP "INPUT HEADING (20 CHARACTERS"
360 DISP "MAX). "
370 INPUT H2$
380 DISP "INPUT DATE (XX/XX/XX). "
390 INPUT H3$
400 PRINT# 1,1 ; N,H2$,H3$
410 GOTO 180
420 REM *****CONSTANT DATA*****
430 CLEAR
440 DISP "INPUT KNOWN MASS OF WGHT#1 (MG). "
450 INPUT A1
460 DISP "INPUT KNOWN MASS OF WGHT#2 (MG). "
470 INPUT B1
480 DISP "INPUT THE STANDARD DEVIATION OF THE BALANCE(MG). "
490 INPUT C1
500 DISP "INPUT THE MASS OF THE SENSITI-"
510 DISP "VITY WGHT. (MG). "
520 INPUT D1
530 DISP "INPUT THE VOLUME OF THE SENSITI-"
540 DISP "VITY WGHT. (CC). "
550 INPUT E1
560 DISP "INPUT TEMPERATURE CORRECTION"
570 DISP "(DEG C). "
580 INPUT F1
590 DISP "INPUT PRESSURE CORRECTION"
600 DISP "(MM HG). "

```

```

610 INPUT G1
620 DISP "INPUT HUMIDITY CORRECTION (%)."
630 INPUT H1
640 PRINT# 1,2 ; A1,B1,C1,D1,E1,F1,G1,H1
650 DISP "INPUT COEF. OF THERM. EXPANSION OF WGHT#1."
660 INPUT A2
670 DISP "INPUT VOLUME OF WGHT#1 (CC)."
680 INPUT B2
690 DISP "INPUT DENSITY OF WGHT#1 (G/CC)."
700 INPUT C2
710 DISP "INPUT MASS OF ADDED WEIGHT TO"
720 DISP "WGHT#1 (MG)."
730 INPUT D2
740 DISP "INPUT COEF. OF THERM. EXPANSION OF THE ADDED WEIGHT."
750 INPUT E2
760 DISP "INPUT VOLUME OF ADDED WEIGHT"
770 DISP "(CC)."
780 INPUT F2
790 PRINT# 1,3 ; A2,B2,C2,D2,E2,F2
800 DISP "INPUT COEF. OF THERM. EXPANSION OF WGHT#2."
810 INPUT A3
820 DISP "INPUT VOLUME OF WGHT#2 (CC)."
830 INPUT B3
840 DISP "INPUT DENSITY OF WGHT#2 (G/CC)."
850 INPUT C3
860 DISP "INPUT MASS OF ADDED WEIGHT TO"
870 DISP "WGHT#2 (MG)."
880 INPUT D3
890 DISP "INPUT COEF. OF THERM. EXPANSION OF ADDED WEIGHT."
900 INPUT E3
910 DISP "INPUT VOLUME OF ADDED WEIGHT"
920 DISP "(CC)."
930 INPUT F3
940 PRINT# 1,4 ; A3,B3,C3,D3,E3,F3
950 DISP "INPUT COEF. OF THERM. EXPANSION OF WGHT#3."
960 INPUT A4
970 DISP "INPUT VOLUME OF WGHT#3 (CC)."
980 INPUT B4
990 DISP "INPUT DENSITY OF WGHT#3 (G/CC)."
1000 INPUT C4
1010 DISP "INPUT MASS OF ADDED WEIGHT TO"
1020 DISP "WGHT#3 (MG)."
1030 INPUT D4
1040 DISP "INPUT COEF. OF THERM. EXPANSION OF ADDED WEIGHT."
1050 INPUT E4
1060 DISP "INPUT VOLUME OF ADDED WEIGHT"
1070 DISP "(CC)."
1080 INPUT F4
1090 PRINT# 1,5 ; A4,B4,C4,D4,E4,F4
1100 DISP "INPUT COEF. OF THERM. EXPANSION OF WGHT#4."
1110 INPUT A5
1120 DISP "INPUT VOLUME OF WGHT#4 (CC)."
1130 INPUT B5
1140 DISP "INPUT DENSITY OF WGHT#4 (G/CC)."
1150 INPUT C5
1160 DISP "INPUT MASS OF ADDED WEIGHT TO"
1170 DISP "WGHT#4 (MG)."
1180 INPUT D5
1190 DISP "INPUT COEF. OF THERM. EXPANSION OF ADDED WEIGHT."
1200 INPUT E5

```

```
1210 DISP "INPUT VOLUME OF ADDED WEIGHT"  
1220 DISP "(CC)."  
1230 INPUT F5  
1240 PRINT# 1,6 ; A5,B5,C5,D5,E5,F5  
1250 GOTO 180  
1260 REM ***** EXIT *****  
1270 ASSIGN# 1 TO *  
1280 SUBEND
```

```

100 SUB "MEASUR" (X7)
110 REM
120 REM *RECORDS BALANCE
130 REM OBSERVATIONS FROM
140 REM THE KEYBOARD.
150 REM
160 REM *RECORDS VOLTMETER
170 REM READINGS OF TEMP,
180 REM PRESS, AND REL HUMD.
190 REM
200 REM *WRITTEN BY J TAYLOR
210 REM FOR OWM 3/85.
220 REM
230 OPTION BASE 1
240 DIM T(6,5),P(6,5),H(6,5),O(6,4),T2(6),P2(6),H2(6)
250 DIM M$(5),X$(2),Z$(13),Y$(2),W$(3)
260 ON X7 GOTO 270,290,310
270 ASSIGN# 1 TO "A MAIN"
280 GOTO 320
290 ASSIGN# 1 TO "B MAIN"
300 GOTO 320
310 ASSIGN# 1 TO "C MAIN"
320 PRINTER IS 2
330 Z$=" INPUT OBSER "
340 W$="+ "
350 ON KEY# 1,"M1:M2" GOTO 490
360 ON KEY# 2,"M1:M3" GOTO 510
370 ON KEY# 3,"M1:M4" GOTO 530
380 ON KEY# 5,"M2:M3" GOTO 550
390 ON KEY# 6,"M2:M4" GOTO 570
400 ON KEY# 7,"M3:M4" GOTO 590
410 ON KEY# 8,"EXIT" GOTO 1080
420 CLEAR @ KEY LABEL
430 DISP USING "4/"
440 DISP "      DEPRESS ANY OF THE"
450 DISP "      LABELED KEYS BELOW TO"
460 DISP "      OBTAIN THE DESIRED"
470 DISP "      MEASUREMENT FUNCTION."
480 GOTO 480
490 M$="M1:M2" @ I=1 @ X$="M1" @ Y$="M2"
500 GOTO 600
510 M$="M1:M3" @ I=2 @ X$="M1" @ Y$="M3"
520 GOTO 600
530 M$="M1:M4" @ I=3 @ X$="M1" @ Y$="M4"
540 GOTO 600
550 M$="M2:M3" @ I=4 @ X$="M2" @ Y$="M3"
560 GOTO 600
570 M$="M2:M4" @ I=5 @ X$="M2" @ Y$="M4"
580 GOTO 600
590 M$="M3:M4" @ I=6 @ X$="M3" @ Y$="M4"
600 CLEAR @ KEY LABEL
610 DISP HGL$(M$) @ PRINT HGL$(M$)
620 PRINT USING "/"
630 Q=1 @ GOSUB 930
640 DISP Z$;X$
650 INPUT O(I,1)
660 DISP Z$;Y$
670 INPUT O(I,2)
680 DISP Z$;Y$;W$
690 INPUT O(I,3)

```

```

700 DISP Z$;X$;W$
710 INPUT O(I,4)
720 Q=2 @ GOSUB 930
730 T1=0 @ P1=0 @ H1=0
740 FOR Q=1 TO 2
750 T1=T1+T(I,Q)
760 P1=P1+P(I,Q)
770 H1=H1+H(I,Q)
780 NEXT Q
790 T1=T1/2 @ P1=P1/2 @ H1=H1/2
800 T2(I)=84.02531*T1-.101
810 P2(I)=(30*P1+800)/1.333224
820 H2(I)=H1/.02-2.5
830 FOR J=1 TO 4
840 PRINT O(I,J)
850 NEXT J
860 PRINT# 1,I+6 ; T2(I),P2(I),H2(I),O(I,1),O(I,2),O(I,3),O(I,4)
870 PRINT USING "/"
880 PRINT USING "4A,2D.3D" ; " T=",T2(I)
890 PRINT USING "4A,3D.3D" ; " P=",P2(I)
900 PRINT USING "4A,2D.2D" ; " H=",H2(I)
910 PRINT USING "2/"
920 GOTO 420
930 REM *** READ AMBIENT ***
940 OUTPUT 724 ; "F1RAN5T2"
950 OUTPUT 710 ; "C1-1+2-3-4-5-6-7-8"
960 GOSUB 1050
970 ENTER 724 ; T(I,Q)
980 OUTPUT 710 ; "C1-1-2+3-4-5-6-7-8"
990 GOSUB 1050
1000 ENTER 724 ; P(I,Q)
1010 OUTPUT 710 ; "C1-1-2-3+4-5-6-7-8"
1020 GOSUB 1050
1030 ENTER 724 ; H(I,Q)
1040 RETURN
1050 WAIT 2000
1060 TRIGGER 724
1070 RETURN
1080 PRINTER IS 701,80
1090 ASSIGN# 1 TO *
1100 SUBEND

```

```

100 SUB "EDITD" (X7)
110 REM
120 REM *EDITS MAIN
130 REM DATA FILES.
140 REM
150 REM *WRITTEN BY J TAYLOR
160 REM FOR OWM 3/85.
170 REM
180 OPTION BASE 1
190 ON X7 GOTO 200,220,240
200 ASSIGN# 1 TO "A MAIN"
210 GOTO 260
220 ASSIGN# 1 TO "B MAIN"
230 GOTO 260
240 ASSIGN# 1 TO "C MAIN"
250 GOTO 260
260 DIM M(8),W(4,6),O(6,4),H2$(20),H3$(8),A(7,3)
270 INTEGER N,L,I,J,K
280 CLEAR
290 DISP USING "2/"
300 DISP "WHICH DATA FILE NEEDS EDITNG"
310 DISP "(1 THRU 5)?"
320 DISP USING "2/"
330 DISP " 1) HEADING/DATE"
340 DISP " 2) CONSTANT DATA"
350 DISP " 3) AMBIENT DATA"
360 DISP " 4) OBSERVATIONS"
370 DISP " 5) EXIT"
380 DISP USING "/".
390 INPUT K
400 IF K<1 OR K>5 THEN GOTO 280
410 ON K GOTO 430,730,1570,1980,2430
420 REM *****HEADING/DATE*****
430 CLEAR
440 READ# 1,1 ; N,H2$,H3$
450 DISP USING "4/"
460 DISP "WHICH PARAMETER NEEDS EDITNG"
470 DISP "(1 THRU 3)?"
480 DISP USING "2/"
490 DISP " 1) HEADING: "
500 DISP " ";H2$
510 DISP
520 DISP " 2) DATE: ";H3$
530 DISP
540 DISP " 3) RETURN"
550 DISP USING "2/"
560 INPUT K
570 IF K<1 OR K>3 THEN GOTO 430
580 ON K GOTO 590,640,280
590 DISP
600 DISP H2$
610 DISP "INPUT NEW HEADING."
620 INPUT H2$
630 GOTO 680
640 DISP
650 DISP H3$
660 DISP "INPUT NEW DATE."
670 INPUT H3$
680 N=2
690 PRINT# 1,1 ; N,H2$,H3$

```

```

700 READ# 1,1 ; N,H2$,H3$
710 GOTO 430
720 REM *****CONSTANT DATA*****
730 CLEAR @ CRT OFF
740 READ# 1,2 ; M(1),M(2),M(3),M(4),M(5),M(6),M(7),M(8)
750 FOR I=1 TO 4
760 READ# 1,I+2 ; W(I,1),W(I,2),W(I,3),W(I,4),W(I,5),W(I,6)
770 NEXT I
780 CRT ON
790 CLEAR
800 DISP USING "/"
810 DISP "WHICH GROUPING NEEDS EDITING"
820 DISP "(1 THRU 6)?"
830 DISP USING "/"
840 DISP " 1) STANDARDS/SENSITIVITY /"
850 DISP "      INSTR. CORREC."
860 DISP " 2) WEIGHT #1"
870 DISP " 3) WEIGHT #2"
880 DISP " 4) WEIGHT #3"
890 DISP " 5) WEIGHT #4"
900 DISP " 6) RETURN"
910 DISP USING "/"
920 INPUT K
930 IF K<1 OR K>6 THEN GOTO 790
940 ON K GOTO 950,1240,1240,1240,1240,280
950 CLEAR
960 DISP USING "3/"
970 DISP " WHICH PARAMETER NEEDS EDITNG"
980 DISP "(1 THRU 9)?"
990 DISP USING "/"
1000 DISP USING 1010 ; "1) MASS WT#1 (MG): ",M(1)
1010 IMAGE 19A,K
1020 DISP USING 1030 ; "2) MASS WT#2 (MG): ",M(2)
1030 IMAGE 19A,K
1040 DISP "3) BALANCE S.D. (MG): ";M(3)
1050 DISP "4) MASS SEN WT (MG):";M(4)
1060 DISP "5) VOL SEN WT (MG): ";M(5)
1070 DISP "6) TEMP CORREC (DEG C): ";M(6)
1080 DISP "7) PRESS CORREC (MM HG): ";M(7)
1090 DISP "8) HUM CORREC (%): ";M(8)
1100 DISP "9) RETURN"
1110 DISP USING "/"
1120 INPUT K
1130 IF K<1 OR K>9 THEN GOTO 950
1140 IF K=9 THEN GOTO 790
1150 DISP
1160 DISP M(K)
1170 DISP "INPUT NEW VALUE."
1180 INPUT M(K)
1190 CRT OFF
1200 PRINT# 1,2 ; M(1),M(2),M(3),M(4),M(5),M(6),M(7),M(8)
1210 READ# 1,2 ; M(1),M(2),M(3),M(4),M(5),M(6),M(7),M(8)
1220 CRT ON @ CLEAR
1230 GOTO 950
1240 K=K-1
1250 CLEAR
1260 DISP USING "3/"
1270 DISP "WHICH PARAMETER NEEDS EDITING"
1280 DISP "(1 THRU 7)?"
1290 DISP USING "2/"

```

```

1300 DISP USING 1310 ; "1) THERM COEF WT#",K,":",W(K,1)
1310 IMAGE X,17A,K,A,X,K
1320 DISP USING 1330 ; "2) VOL OF WT#",K,":",W(K,2)
1330 IMAGE X,13A,K,A,X,K
1340 DISP USING 1350 ; "3) DEN OF WT#",K,":",W(K,3)
1350 IMAGE X,13A,K,A,X,K
1360 DISP USING 1370 ; "4) MASS OF ADD WT:",W(K,4)
1370 IMAGE X,18A,X,K
1380 DISP USING 1390 ; "5) THERM COEF ADD WT:",W(K,5)
1390 IMAGE X,21A,X,K
1400 DISP USING 1410 ; "6) VOL OF ADD WT:",W(K,6)
1410 IMAGE X,17A,X,K
1420 DISP " 7) RETURN"
1430 DISP USING "/"
1440 INPUT L
1450 IF L<1 OR L>7 THEN GOTO 1260
1460 IF L=7 THEN GOTO 790
1470 DISP
1480 DISP W(K,L)
1490 DISP "INPUT NEW VALUE"
1500 INPUT W(K,L)
1510 CRT OFF @ CLEAR
1520 PRINT# 1,K+2 ; W(K,1),W(K,2),W(K,3),W(K,4),W(K,5),W(K,6)
1530 READ# 1,K+2 ; W(K,1),W(K,2),W(K,3),W(K,4),W(K,5),W(K,6)
1540 CRT ON
1550 GOTO 1250
1560 REM *****AMBIENT DATA*****
1570 CLEAR @ CRT OFF
1580 FOR I=1 TO 6
1590 READ# 1,I+6 ; A(I,1),A(I,2),A(I,3),O(I,1),O(I,2),O(I,3),O(I,4)
1600 NEXT I
1610 CRT ON
1620 CLEAR
1630 DISP USING "2/"
1640 DISP "WHICH DATA POINT NEEDS EDITING"
1650 DISP "(1 THRU 7)?"
1660 DISP USING "/"
1670 T$=HGL$("T")
1680 P$=HGL$("P")
1690 R$=HGL$("RH")
1700 DISP USING 1710 ; T$,P$,R$
1710 IMAGE 8X,A,7X,A,6X,2A,2/
1720 FOR I=1 TO 6
1730 DISP USING 1740 ; I,")",A(I,1),A(I,2),A(I,3)
1740 IMAGE 2X,K,A,2X,DZ.2D,2X,2DZ.3D,2X,DZ.D
1750 NEXT I
1760 DISP USING "/"
1770 DISP " 7) RETURN"
1780 INPUT K
1790 IF K<1 OR K>7 THEN GOTO 1620
1800 IF K=7 THEN GOTO 280
1810 DISP
1820 DISP USING 1830 ; K,")",A(K,1),A(K,2),A(K,3)
1830 IMAGE X,K,A,X,K,4X,K,4X,K
1840 DISP "WHICH VALUE NEEDS EDITING"
1850 DISP "(1,2,3)?"
1860 INPUT L
1870 IF L<1 OR L>3 THEN GOTO 1810
1880 DISP
1890 DISP A(K,L)

```

```

1900 DISP "INPUT NEW VALUE"
1910 INPUT A(K,L)
1920 CRT OFF
1930 PRINT# 1,K+6 ; A(K,1),A(K,2),A(K,3),O(K,1),O(K,2),O(K,3),O(K,4)
1940 READ# 1,K+6 ; A(K,1),A(K,2),A(K,3),O(K,1),O(K,2),O(K,3),O(K,4)
1950 CRT ON @ CLEAR
1960 GOTO 1630
1970 REM *****OBSERVATIONS*****
1980 CLEAR @ CRT OFF
1990 FOR I=1 TO 6
2000 READ# 1,I+6 ; A(I,1),A(I,2),A(I,3),O(I,1),O(I,2),O(I,3),O(I,4)
2010 NEXT I
2020 CRT ON
2030 CLEAR
2040 DISP USING "/"
2050 DISP "WHICH COMPARISON NEEDS EDITING"
2060 DISP "(1 THRU 7)?"
2070 DISP USING "/"
2080 DISP " 1) WGHT#1 - WGHT#2"
2090 DISP " 2) WGHT#1 - WGHT#3"
2100 DISP " 3) WGHT#1 - WGHT#4"
2110 DISP " 4) WGHT#2 - WGHT#3"
2120 DISP " 5) WGHT#2 - WGHT#4"
2130 DISP " 6) WGHT#3 - WGHT#4"
2140 DISP " 7) RETURN"
2150 DISP USING "/"
2160 INPUT K
2170 IF K<1 OR K>7 THEN GOTO 2030
2180 IF K=7 THEN GOTO 280
2190 CLEAR
2200 DISP USING "4/"
2210 DISP "WHICH PARAMETER NEEDS EDITING"
2220 DISP "(1 THRU 5)?"
2230 DISP USING "2/"
2240 FOR I=1 TO 4
2250 DISP USING 2260 ; I,")",O(K,I)
2260 IMAGE 6X,K,A,X,K
2270 NEXT I
2280 DISP " 5) RETURN"
2290 DISP USING "2/"
2300 INPUT L
2310 IF L<1 OR L>5 THEN GOTO 2190
2320 IF L=5 THEN GOTO 2030
2330 DISP
2340 DISP O(K,L)
2350 DISP "INPUT NEW VALUE"
2360 INPUT O(K,L)
2370 CRT OFF @ CLEAR
2380 PRINT# 1,K+6 ; A(K,1),A(K,2),A(K,3),O(K,1),O(K,2),O(K,3),O(K,4)
2390 READ# 1,K+6 ; A(K,1),A(K,2),A(K,3),O(K,1),O(K,2),O(K,3),O(K,4)
2400 CRT ON
2410 GOTO 2190
2420 REM ***** EXIT *****
2430 ASSIGN# 1 TO *
2440 SUBEND

```

```

100 SUB "ANALY" (X7)
110 REM
120 REM *COMPUTES TRUE MASS
130 REM USING FOUR 1'S SERIES
140 REM
150 REM *DERIVED FROM FTM77
160 REM PROGRAM BY C L CARROLL
170 REM AND R S DAVIS.
180 REM
190 REM *WRITTEN BY J TAYLOR
200 REM FOR OWM 3/85.
210 REM
220 OPTION BASE 1
230 REAL A(7,3),B(6),C(4,2),D(4),E(4,4),F(4,6),G(6),H(6),O(6),P(4)
240 REAL Q(6),R(6),S(6),T(6),V(4,2),W(4),Y(4),Z(4),A5,B5,B1,T1
250 REAL V2,F5,M5,M1,M2,M3,M4,O1(5),Q3,D1,S1,S2,S3,V1(6),X1,Y1,H1
260 REAL T2,R2,T4,Q1,Q2,Z1,T3,R3
270 INTEGER L1(7,4),L2(7,4),M(6,6),N,L3(7,4),L4(7,4),L(7,4)
280 DIM H2$(20),H3$(8)
290 FOR I=1 TO 4
300 FOR J=1 TO 7
310 READ L1(J,I)
320 NEXT J
330 NEXT I
340 FOR I=1 TO 4
350 FOR J=1 TO 7
360 READ L2(J,I)
370 NEXT J
380 NEXT I
390 FOR I=1 TO 4
400 FOR J=1 TO 7
410 READ L3(J,I)
420 NEXT J
430 NEXT I
440 FOR I=1 TO 4
450 FOR J=1 TO 7
460 READ L4(J,I)
470 NEXT J
480 NEXT I
490 DATA 0,0,0,0,0,0,4,-2,-1,-1,1,1,0,4
500 DATA -1,-2,-1,-1,0,1,4,-1,-1,-2,0,-1,-1,4
510 DATA 2,1,1,-1,-1,0,4,-2,-1,-1,1,1,0,4
520 DATA 0,-3,-1,-3,-1,2,4,0,-1,-3,-1,-3,-2,4
530 DATA 3,3,2,0,-1,-1,4,-3,0,-1,3,2,-1,4
540 DATA 0,-3,-1,-3,-1,2,4,0,0,-4,0,-4,-4,4
550 DATA 1,1,1,0,0,0,1,-1,0,0,1,1,0,1
560 DATA 0,-1,0,-1,0,1,1,0,0,-1,0,-1,-1,1
570 FOR I=1 TO 6
580 FOR J=1 TO 6
590 READ M(J,I)
600 NEXT J
610 NEXT I
620 DATA 4,-2,-2,2,2,0,-2,4,-2,-2,0,2,-2,-2,4,0,-2,-2
630 DATA 2,-2,0,4,-2,2,2,0,-2,-2,4,-2,0,2,-2,2,-2,4
640 FOR I=1 TO 4
650 FOR J=1 TO 4
660 READ E(J,I)
670 NEXT J
680 NEXT I
690 DATA 0,.5,.5,.5,.125,.125,.375,.375,.166666666667,.166666666667,.166666666667

```

```

7
700 DATA .3333333333333333,.1875,.1875,.1875,.1875
710 F5=3.78
720 A5=3
730 P(1)=4
740 P(2)=8
750 P(3)=12
760 P(4)=4
770 ON X7 GOTO 780,800,820
780 ASSIGN# 1 TO "A MAIN"
790 GOTO 840
800 ASSIGN# 1 TO "B MAIN"
810 GOTO 840
820 ASSIGN# 1 TO "C MAIN"
830 GOTO 840
840 READ# 1,1 ; N,H2$,H3$
850 READ# 1,2 ; M1,M2,B5,S3,S2,T1,B1,H1
860 M3=0
870 M4=0
880 FOR I=1 TO 4
890 READ# 1,I+2 ; C(I,1),V(I,1),D(I),W(I),C(I,2),V(I,2)
900 NEXT I
910 T2=0
920 R2=0
930 PRINT USING 940 ; H2$,H3$
940 IMAGE 30X,20A,20X,8A,3/
950 PRINT USING 960 ; "OBS 1","OBS 2","OBS 3","OBS 4","MG/DV","DRIFT"
960 IMAGE 5X,K,5X,K,5X,K,5X,K,7X,K,6X,K,/
970 FOR I=1 TO 6
980 READ# 1,I+6 ; T(I),B(I),H(I),O1(1),O1(2),O1(3),O1(4)
990 T(I)=T(I)+T1
1000 T2=T2+T(I)
1010 B(I)=B(I)+B1
1020 H(I)=H(I)+H1
1030 V1(I)=0
1040 R(I)=0
1050 X2=V1(I)
1060 X3=R(I)
1070 X4=T(I)
1080 X5=B(I)
1090 X6=H(I)
1100 CALL "SUBRHO" ( X4,X5,X6,X2,X3 )
1110 V1(I)=X2
1120 R(I)=X3
1130 R2=R2+R(I)
1140 T(I)=T(I)-20
1150 S1=O1(1)-O1(4)+3*(O1(3)-O1(2))
1160 S(I)=2*(S3-R(I)*S2)/S1
1170 A(I,1)=(O1(1)-O1(2)-O1(3)+O1(4))*S(I)/2
1180 G(I)=(-O1(1)+O1(2)-O1(3)+O1(4))/2
1190 G(I)=G(I)*S(I)
1200 PRINT USING 1210 ; O1(1),O1(2),O1(3),O1(4),S(I),G(I)
1210 IMAGE X,4(M4DZ.3D),2X,M2DZ.5D,2X,M2DZ.3D
1220 NEXT I
1230 Q(1)=W(1)-W(2)
1240 Q(2)=W(1)-W(3)
1250 Q(3)=W(1)-W(4)
1260 Q(4)=W(2)-W(3)
1270 Q(5)=W(2)-W(4)
1280 Q(6)=W(3)-W(4)

```

```

1290 FOR I=1 TO 2
1300 F(I,1)=V(1,I)*(1+C(1,I)*T(1))-V(2,I)*(1+C(2,I)*T(1))
1310 F(I,2)=V(1,I)*(1+C(1,I)*T(2))-V(3,I)*(1+C(3,I)*T(2))
1320 F(I,3)=V(1,I)*(1+C(1,I)*T(3))-V(4,I)*(1+C(4,I)*T(3))
1330 F(I,4)=V(2,I)*(1+C(2,I)*T(4))-V(3,I)*(1+C(3,I)*T(4))
1340 F(I,5)=V(2,I)*(1+C(2,I)*T(5))-V(4,I)*(1+C(4,I)*T(5))
1350 F(I,6)=V(3,I)*(1+C(3,I)*T(6))-V(4,I)*(1+C(4,I)*T(6))
1360 NEXT I
1370 FOR I=1 TO 6
1380 A(I,2)=A(I,1)+R(I)*(F(1,I)+F(2,I))-Q(I)
1390 NEXT I
1400 A(7,2)=M1
1410 IF N=2 THEN A(7,2)=M1+M2
1420 A(1,3)=M1-(A(1,2)+R(1)*V(2,1)*(1+C(2,1)*T(1)))
1430 A(2,3)=M1-(A(2,2)+R(2)*V(3,1)*(1+C(3,1)*T(2)))
1440 A(3,3)=M1-(A(3,2)+R(3)*V(4,1)*(1+C(4,1)*T(3)))
1450 A(4,3)=M2-(A(4,2)+R(4)*V(3,1)*(1+C(3,1)*T(4)))
1460 A(5,3)=M2-(A(5,2)+R(5)*V(4,1)*(1+C(4,1)*T(5)))
1470 A(6,3)=0
1480 PRINT USING 1490
1490 IMAGE //
1500 PRINT USING 1510 ; "TEMP","BARO","HUM","VAP","RHO","DV","RHO*DV"
1510 IMAGE 6X,K,5X,K,5X,K,6X,K,6X,K,8X,K,7X,K,/
1520 FOR I=1 TO 6
1530 T4=T(I)+20
1540 V2=F(1,I)+F(2,I)
1550 D1=R(I)*V2
1560 PRINT USING 1570 ; T4,B(I),H(I),V1(I),R(I),V2,D1
1570 IMAGE 3X,M2DZ.2D,X,M4DZ.2D,M4DZ.1D,M3DZ.2D,X,M2DZ.4D,M4DZ.4D,X,M4DZ.3D
1580 NEXT I
1590 PRINT USING 1600
1600 IMAGE 2/
1610 PRINT USING 1620 ; "BAL DIFF","MASS IN AIR","A","DELTA"
1620 IMAGE 6X,K,5X,K,13X,K,12X,K,/
1630 Q3=0
1640 FOR I=1 TO 6
1650 X1=0
1660 FOR J=1 TO 6
1670 X1=X1+M(J,I)*A(J,2)
1680 NEXT J
1690 O(I)=X1/8
1700 Q3=Q3+O(I)^2
1710 PRINT USING 1720 ; A(I,1),A(I,3),A(I,2),O(I)
1720 IMAGE 5X,M2DZ.3D,2X,M9DZ.3D,X,M9DZ.3D,4X,M5DZ.3D
1730 NEXT I
1740 PRINT USING 1750
1750 IMAGE 2/
1760 PRINT USING 1770 ; "DENSITY","TM","SD","V AT 20"
1770 IMAGE 7X,K,9X,K,17X,K,10X,K,/
1780 FOR I=1 TO 4
1790 Y1=0
1800 FOR J=1 TO 7
1810 IF N=1 THEN L(J,I)=L1(J,I)
1820 IF N=2 THEN L(J,I)=L2(J,I)
1830 Y1=Y1+L(J,I)*A(J,2)
1840 NEXT J
1850 Y(I)=Y1/P(N)
1860 V(I,1)=Y(I)*.001/D(I)
1870 IF E(I,N)=0 THEN 1910
1880 Z(I)=SQRT(B5^2*E(I,N))

```

```

1900 GOTO 1920
1910 Z(I)=0
1920 PRINT USING 1930 ; D(I),Y(I),Z(I),V(I,1)
1930 IMAGE 4X,M2DZ.5D,X,M9DZ.3D,,X,M9DZ.3D,3X,M5DZ.5D
1940 NEXT I
1950 Q1=SQR(Q3/A5)
1960 Q2=Q3/A5/B5^2
1970 PRINT USING 1980
1980 IMAGE 2/
1990 PRINT USING 2000 ; "ACC SD","OBS SD","F RATIO","F TEST","DF"
2000 IMAGE 7X,K,8X,K,6X,K,6X,K,9X,K,/
2010 PRINT USING 2020 ; B5,Q1,Q2,F5,A5
2020 IMAGE 2(M7DZ.3D,X),3X,M3DZ.3D,3X,M4DZ.2D,10X,D
2030 IF ABS(Q2)<=3.78 THEN GOTO 2050
2040 PRINT "***** WARNING ***** F RATIO IS NOT IN CONTROL."
2050 PRINT USING 1980
2060 IF N=1 THEN M5=M2
2070 IF N=2 THEN M5=M1-M2
2080 IF N=1 THEN Y1=Y(2)
2090 IF N=2 THEN Y1=Y(1)-Y(2)
2095 IF N=2 THEN Z(2)=SQR(B5^2*.5)
2100 Z1=(Y1-M5)/Z(2)
2110 T3=T2/6
2120 R3=R2/6
2130 PRINT USING 2140 ; "ACC DIFF","OBS DIFF","T TEST","MEAN TEMP","MEAN RHO"
2140 IMAGE 7X,K,6X,K,7X,K,4X,K,4X,K,/
2150 PRINT USING 2160 ; M5,Y1,Z1,T3,R3
2160 IMAGE 3X,2(M6DZ.3D,2X),X,M4DZ.3D,3X,M3DZ.2D,3X,M3DZ.4D
2170 IF ABS(Z1)<=3 THEN GOTO 2190
2180 PRINT "***** WARNING ***** T TEST IS NOT IN CONTROL."
2190 ASSIGN# 1 TO *
2200 SUBEND

```

```

100 SUB "SUBRHO" (T1,P1,R1,V,R2)
110 REM
120 REM *CALCULATES BOUYANCY
130 REM CORRECTION.
140 REM
150 REM *DERIVED FROM FTN77
160 REM PROGRAM BY R S DAVIS
170 REM
180 REM *WRITTEN BY J TAYLOR
190 REM FOR OWM 3/85
200 REM
210 P2=.043
220 A3=(28.9635+12.011*(P2/100-.0004))*0.001
230 R=8.31441
240 C2=A3/R
250 T=273.15+T1
260 P=P1*13.5951*9.80665
270 A=.000012811805
280 C=34.04926034
290 B=-.019509874
300 D=-6353.6311
310 P3=EXP(A*T*T+B*T+C+D/T)
320 A4=1.00062
330 B2=.0000000314
340 G=.000000056
350 F=A4+B2*P+G*T1*T1
360 H=R1/100
370 X=H*F*P3/P
380 V=X*P1
390 A0=.00000162419
400 A2=1.088E-10
410 B0=.0000005757
420 C0=.00019297
430 D0=1.73E-11
440 A1=-.000000028969
450 B1=-.00000002589
460 C1=-.000002285
470 E=-.00000001034
480 Z=1-P/T*(A0+A1*T1+A2*T1*T1+(B0+B1*T1)*X+(C0+C1*T1)*X*X)+P^2/T^2*(D0+E*X*X)
490 R2=C2*P/(T*Z)*(1-.378*X)
500 SUBEND

```

```

00 SUB "PRINTM" (X7)
10 REM
20 REM *PRINTS CONTENTS OF THE
30 REM MAIN DATA FILE.
40 REM
50 REM *WRITTEN BY J TAYLOR
60 REM FOR OWM 3/85.
70 REM
80 ON X7 GOTO 190,210,230
90 ASSIGN# 1 TO "A MAIN"
00 GOTO 250
10 ASSIGN# 1 TO "B MAIN"
20 GOTO 250
30 ASSIGN# 1 TO "C MAIN"
40 GOTO 250
50 DIM H2$(20),H3$(8),T$(24)
60 READ# 1,1 ; N,H2$,H3$
70 PRINT USING "6/"
80 T$="**** MAIN DATA FILE ****"
90 PRINT USING "25X,24A,2/" ; T$
00 PRINT USING 330 ; N,H2$,H3$
10 PRINT USING 320
20 IMAGE /
30 IMAGE D,5X,20A,30X,8A
40 READ# 1,2 ; A,B,C,D,E,F,G,H
50 PRINT USING 360 ; A,B,C,D,E,F,G,H
60 IMAGE 2(6DZ.3D,2X),Z.3D,X,2DZ.5D,X,DZ.5D,X,3(DZ.3D,X)
70 PRINT USING 320
80 FOR I=1 TO 4
90 READ# 1,I+2 ; A,B,C,D,E,F
00 PRINT USING 410 ; A,B,C,D,E,F
10 IMAGE Z.6D,3X,2DZ.4D,3X,2DZ.5D,3X,5DZ.6D,3X,Z.6D,3X,DZ.5D
20 NEXT I
30 PRINT USING 320
40 FOR I=1 TO 6
50 READ# 1,I+6 ; A,B,C,D,E,F,G
60 PRINT USING 470 ; A,B,C,D,E,F,G
70 IMAGE 8(M2DZ.3D,2X)
80 NEXT I
90 ASSIGN# 1 TO *
00 SUBEND

```

```

100 REM PROGRAM "TRANSD"
110 REM
120 REM *TRANSFERS DATA FILE
130 REM FROM AN OLD TAPE TO
140 REM A NEW TAPE CARTRIDGE.
150 REM
160 REM *WRITTEN BY J TAYLOR
170 REM FOR OWM 3/85.
180 REM
190 OPTION BASE 1
200 DISP "INPUT NAME OF THE DATA FILE."
210 INPUT A$
220 ASSIGN# 1 TO A$
230 DIM H2$(20),H3$(8),A(11),B(11),C(11),D(11),E(11),F(11),G(11),H(11)
240 INTEGER N
250 READ# 1,1 ; N,H2$,H3$
260 READ# 1,2 ; A(1),B(1),C(1),D(1),E(1),F(1),G(1),H(1)
270 FOR I=1 TO 4
280 READ# 1,I+2 ; A(I+1),B(I+1),C(I+1),D(I+1),E(I+1),F(I+1)
290 NEXT I
300 FOR I=1 TO 6
310 READ# 1,I+6 ; A(I+5),B(I+5),C(I+5),D(I+5),E(I+5),F(I+5),G(I+5)
320 NEXT I
330 ASSIGN# 1 TO *
340 CLEAR
350 DISP
360 DISP
370 DISP
380 DISP
390 DISP " 1)REMOVE OLD TAPE."
400 DISP
410 DISP " 2)INSERT NEW TAPE."
420 DISP
430 DISP " 3)PURGE DATA FILE IN NEW"
440 DISP " TAPE,IF NECESSARY."
450 DISP
460 DISP " 4)RESUME BY PRESSING CONT."
470 PAUSE
480 CREATE A$,12,80
490 ASSIGN# 1 TO A$
500 PRINT# 1,1 ; N,H2$,H3$
510 PRINT# 1,2 ; A(1),B(1),C(1),D(1),E(1),F(1),G(1),H(1)
520 FOR I=1 TO 4
530 PRINT# 1,I+2 ; A(I+1),B(I+1),C(I+1),D(I+1),E(I+1),F(I+1)
540 NEXT I
550 FOR I=1 TO 6
560 PRINT# 1,I+6 ; A(I+5),B(I+5),C(I+5),D(I+5),E(I+5),F(I+5),G(I+5)
570 NEXT I
580 ASSIGN# 1 TO *
590 CLEAR
600 CAT
610 END

```

Appendix III.

Once a laboratory has recorded a complete set of data onto the tape cartridge an analysis of this information can begin and the 80-column printer will automatically generate a single page of tabulated results followed by a one-half page printout of the information stored in the data file. A sample of this output is shown in table 1 using data stored in one of the three data files, table 2. The assigned true mass values, TM, of each kilogram are located in the fourth row, second column of the printout. Where the first two of four values in this group correspond to assigned true mass values of the laboratory standards, and the third and fourth values in this group are the true mass values assigned to the "unknown" artifacts (A1-S2, R2-H2, or A1-TE).

The following is a brief description of each quality found in this report:

OBS (x) - Refers to each of the four balance observations of the double substitution weighing. Observations are usually read in milligrams; however, the NIST balance with its servo generates values in volts. The nature of substitution weighing with a sensitivity weight is that the observations may be of an arbitrary scale but in the data reduction the proper unit is established.

MG/DV - The value assigned to the balance optical scale divisions from weighing a sensitivity weight.

DRIFT - Linear shift (mg) between balance observations.

TEMP - Average temperature ($^{\circ}\text{C}$) obtained from automatic readings taken before and after a double substitution weighing.

BARO - Average barometric pressure (mm Hg) obtained from automatic readings taken before and after a double substitution weighing.

HUM - Average relative humidity (%) obtained from automatic readings taken before and after a double substitution weighing.

VAP - Vapor pressure (mm Hg) of water.

RHO - Computed air density (mg/cm^3) derived by subprogram SUBRHO. This calculation is based on the work of Jones [5].

DV - Difference in volume (cm^3) between kilograms compared in the double substitution weighing.

RHO*DV - Mass (mg) correction that arises from the difference in buoyant forces acting upon the two kilograms.

BAL DIFF - Difference in the observed readings (mg) between two kilograms. This mass difference includes any tare weight that may be added to the kilograms, but corrections due to the effects of air buoyancy are not applied to this group of weights.

MASS IN AIR - Assigned mass (mg) without buoyancy corrections, of the "unknown" kilogram or the check standard with respect to the laboratory standard used in the double substitution weighing.

A - Difference in the observed readings (mg) between two kilograms. In this calculation the mass of all added tare weights have been removed and the appropriate air buoyancy correction applied to each artifact.

DELTA - Residuals (mg) obtained when fitting the data using the method of least squares.

DENSITY - Density (g/cm³) of each kilogram based on hydrostatic measurement.

TM - The mass (mg) of each kilogram.

SD - The measurement estimate of standard deviation (mg) for each kilogram.

V AT 20 - Volume (cm³) of each kilogram at 20 degrees Celsius. This group of values is derived from the computed true mass values and densities of the kilograms.

ACC SD - Accepted standard deviation (mg) of the balance.

OBS SD - The observed estimate of standard deviation (mg) of the balance as determined for each weighing series by the least squares adjustments.

F RATIO - Examines the precision of the weighing series using the equation

$$F \text{ ratio} = \frac{(\text{OBS SD})^2}{(\text{ACC SD})^2}$$

where,

OBS SD = observed standard deviation of the balance, and
ACC SD = accepted standard deviation of the balance.

F TEST - Has a constant value equal to 3.78 (the critical value) and is compared to the weighing series F ratio. A warning is printed if the absolute value of F ratio is greater than the F test. The value 3.78 is based on the number of degrees of freedom in OBS SD (3) and ACC SD (∞).

DF - Degrees of freedom.

ACC DIFF - The accepted difference in mass (mg) of the two laboratory standards.

OBS DIFF - The measured difference in mass (mg) of the two laboratory standards.

t-TEST - Examines the difference in mass between the two standards using the equations below. However, the equation used here ignores the random uncertainty associated with the standards. In a calibration laboratory, the random uncertainty is small due to averaging over a large pool of data.

$$t\text{-test} = \frac{(\text{ACC DIFF}) - (\text{OBS DIFF})}{S}$$

and

$$S = [(\text{ACC SD})^2/2]^{1/2}$$

where,

ACC DIFF = accepted difference between the two laboratory standards,

OBS DIFF = observed difference between the two laboratory standards,

and

ACC SD = accepted standard deviation of the balance.

MEAN TEMP - Average temperature (°C) of the six double substitutions.

MEAN RHO - Average air density (mg/cm³) of the six double substitutions.

Table 1. Printout of results derived from data file in table 2.

SAMPLE DATA					2/5/86	
OBS 1	OBS 2	OBS 3	OBS 4	MG/DV	DRIFT	
22.600	24.270	43.900	42.440	1.02418	0.108	
22.850	21.080	40.690	42.560	1.02235	0.051	
22.980	21.850	41.460	42.730	1.02339	0.072	
24.660	21.350	40.960	44.300	1.02052	0.015	
24.680	21.900	41.560	44.280	1.01560	-0.030	
21.480	21.940	41.520	41.120	1.02287	0.031	
TEMP	BARO	HUM	VAP	RHO	DV	RHO*DV
23.66	743.43	38.8	8.56	1.1586	-0.2398	-0.278
23.67	743.44	38.8	8.55	1.1587	-232.1703	-269.005
23.67	743.47	38.8	8.56	1.1587	67.3729	78.063
23.67	743.48	38.8	8.56	1.1587	-231.9306	-268.730
23.68	743.47	38.8	8.56	1.1586	67.6128	78.339
23.68	743.48	38.8	8.57	1.1587	299.5435	347.067
BAL DIFF		MASS IN AIR		A	DELTA	
-1.603		999865.689		-1.881	-0.028	
1.861		999862.223		-267.144	0.067	
1.228		999857.860		84.287	-0.038	
3.393		999862.238		-265.337	0.021	
2.793		999857.848		86.128	-0.049	
-0.440		0.000		351.623	0.088	
DENSITY	TM		SD		V AT 20	
7.85030	1000011.691		0.018		127.38516	
7.83556	1000013.543		0.018		127.62503	
2.78252	1000278.901		0.031		359.48669	
16.65790	999927.366		0.031		60.02722	
ACC SD	OBS SD	F RATIO		F TEST	DF	
0.050	0.076	2.306		3.78	3	
ACC DIFF	OBS DIFF	T TEST		MEAN TEMP	MEAN RHO	
-1.828	-1.852	-0.688		23.67	1.1587	

Table 2. Printout of data file used to generate results in table 9.

**** MAIN DATA FILE ****

2	SAMPLE DATA						2/5/86
1000011.703	1000013.531	0.050	20.00570	0.00741	0.000	0.000	0.000
0.000045	127.3852	7.85030		0.000000	0.000000	0.00000	
0.000045	127.6250	7.83556		0.000000	0.000000	0.00000	
0.000069	359.4856	2.78252		0.000000	0.000000	0.00000	
0.000020	60.0271	16.65790		4.995900	0.000069	0.00185	
23.663	743.428	38.847	22.600	24.270	43.900	42.440	
23.666	743.444	38.792	22.850	21.080	40.690	42.560	
23.670	743.467	38.846	22.980	21.850	41.460	42.730	
23.673	743.476	38.847	24.660	21.350	40.960	44.300	
23.676	743.469	38.846	24.680	21.900	41.560	44.280	
23.679	743.484	38.847	21.480	21.940	41.520	41.120	

U.S. DEPT. OF COMM. BIBLIOGRAPHIC DATA SHEET <i>(See instructions)</i>	1. PUBLICATION OR REPORT NO. NISTIR 88-3876	2. Performing Organ. Report No.	3. Publication Date NOVEMBER 1988
4. TITLE AND SUBTITLE <p style="text-align: center;">An Investigation of a User-Operated Mass Calibration Package</p>			
5. AUTHOR(S) Randall M. Schoonover and James E. Taylor			
6. PERFORMING ORGANIZATION <i>(If joint or other than NBS, see instructions)</i> <p style="text-align: center;">NATIONAL BUREAU OF STANDARDS U.S. DEPARTMENT OF COMMERCE GAITHERSBURG, MD 20899</p>			7. Contract/Grant No. 8. Type of Report & Period Covered
9. SPONSORING ORGANIZATION NAME AND COMPLETE ADDRESS <i>(Street, City, State, ZIP)</i>			
10. SUPPLEMENTARY NOTES <input type="checkbox"/> Document describes a computer program; SF-185, FIPS Software Summary, is attached.			
11. ABSTRACT <i>(A 200-word or less factual summary of most significant information. If document includes a significant bibliography or literature survey, mention it here)</i> <p>Reported here are the results of two round-robin mass measurement programs. The first round-robin elucidated the many technical problems that required solutions before successful mass calibrations could be performed in all of the participant laboratories. This report details the technical innovations i.e., thermal conditioning of the kilogram weights, balance servo control, automatic data acquisition, the measurement of some air density parameters, computer software, and presents the results. We believe the results clearly indicate we could, in the future, successfully calibrate mass standards at locations remote from the NBS laboratory while maintaining the rigor necessary for certification.</p>			
12. KEY WORDS <i>(Six to twelve entries; alphabetical order; capitalize only proper names; and separate key words by semicolons)</i> air buoyancy; balance servo; computer-control; mass; software; surface-effects; temperature control; thermal-conditioning; thermal-effects; weighing			
13. AVAILABILITY <input checked="" type="checkbox"/> Unlimited <input type="checkbox"/> For Official Distribution. Do Not Release to NTIS <input type="checkbox"/> Order From Superintendent of Documents, U.S. Government Printing Office, Washington, D.C. 20402. <input checked="" type="checkbox"/> Order From National Technical Information Service (NTIS), Springfield, VA. 22161			14. NO. OF PRINTED PAGES <p style="text-align: center;">64</p> 15. Price <p style="text-align: center;">\$13.95</p>

